


RDNS EQUIPMENT CATALOGUE

EXTERNAL PRESCRIBERS

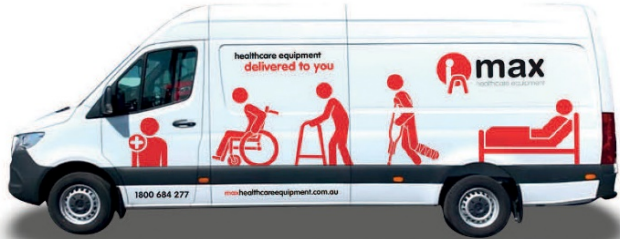

CHSP EQUIP
IN SCOPE LIST

PRODUCT

PRODUCT DESCRIPTION

SPECIFICATIONS

IMAGE

BATHING AND TOILETING

Bath Aids

Yes

Bathboard with handle

Provides a support surface to place across the bath which enables safe and easy showering in a seated position. Helps when transferring in and out of the bath as this can be performed in a seated position. The included handle assists with transfers in and out of the bath as well as safety whilst seated showering.

Adjustable width clamps:
41 cm - 62 cm
Length: 69 cm
Width: 26 cm
Maximum user weight capacity: 150 kg





Yes

Bathboard - Bariatric

Designed to provide a stable seat over the bath.
Constructed of reinforced PVC
Anti - slip rubber pads improves safety and protects the bath surface

Adjustable width clamps:
50 cm - 66 cm
Length: 66 cm
Width: 34 cm
Maximum user weight capacity: 190 kg



Yes	Bath Seat - Swivel	Designed for showering over bath and aiding with transfers into the bath by allowing this occur in a safe seated position.	<p>Maximum user weight: 110kg</p> <p>Adjustable width to fit most baths</p> <p>Moulded plastic seat with drainage holes</p> <p>Contoured back rest</p> <p>Side support arms for added security</p> <p>360° swivel seat</p> <p>Locking lever for set positioning</p>	
Yes	Bath Transfer Bench - With Back	<p>Provides a support surface to place across the bath which enables safe and easy showering in a seated position. Helps when transferring in and out of the bath as this can be performed in a seated position.</p> <p>Plastic seat and backrest provides a non slip surface with holes for water drainage</p> <p>Adjustable backrest and handle to suit left or right facing baths without using tools</p> <p>Built-in dish and hand-held shower slot for user convenience</p>	<p>Height adjustable seat: 43 cm - 56 cm (25mm increments) Seat dimension: Width 66 cm Depth 37 cm Maximum user weight: 136kg</p>	

Yes

Shower Mat / Bath Mat Non Slip

With their raised ridges and unique design, the mats allow you to stand, without slipping, on a wet surface. Use on flat firm wet surface & press down first so the suction caps engage. The mats grip securely and safely due to their high-grade medical rubber content and additional suction cups.

Machine washable at 40°.
Dimension: 555mm x 555mm



Mobile Shower Commodes



No

Mobile Shower Commode,
Transit, 457mm/18inch Wide,
Aluminium

The attendant propelled mobile shower commode is lightweight and maneuverable to allow for easier toileting and showering whilst the user is seated. The chair features a lightweight, silver painted aluminium frame, which provides a rust free, high strength product. Brakes on all four wheels for security, the swing away, detachable footrests offer comfort and support whilst being transported and the armrests can be flipped back to facilitate side transfer. This chair can be used in a shower cubicle or wet room, used as a commode or over a standard toilet bowl. The front of the seat is cut-away to facilitate personal hygiene and a set insert is provided for comfort during transport. Supplied with a square commode pan as standard. The detachable back enables the shower commode chairs to be partially dismantled for more compact storage.

Overall Height: 958 –
1008mm (with castors)
Overall Length: 550mm
(without footrests) / 810mm
(with footrests)
Overall Width: 555mm
Seat Height: 495 – 570mm
Seat Width: 457mm Seat
Depth: 425mm
Aperture Size: 190mm x
360mm
Backrest Width: 462mm
Backrest Height (from
seat): 436mm
Front & Rear Wheel
Diameter: 130mm
Under Seat Clearance: 430
– 505mm
Weight: 12.4kg
Maximum User Weight:
125kg
Warranty: 12 months



No	Mobile Shower Commode, Transit, Stainless Steel, 460mm	<p>Assists with the transportation of the user for showering and toileting. The chair can be wheeled over a standard toilet and used in a large shower alcove. Swingback arms with padded armpads assist side transfers and can be positioned at two different heights Adjustable safety arms for additional security Sliding footrest is textured and weight-bearing for added safety Footrest slides back under frame to assist patient transfers PVC slatted backrest Mounted on 4 x 125mm stainless steel locking castors Durable stainless steel frame Tested and conforms to Australian Standards AS/NZS 3973: 2009I</p>	<p>Seat Width - 460mm Seat Depth - 430mm Seat Height - 595 Armrest Height - 200/240mm Backrest Height - 340mm Overall Width - 540mm Max User Weight - 200kg</p>	
No	Etac Swift Mobile Shower Commode, Transit, Swing Away Footplates, Adjustable Height	<p>Assists with the transportation of the user for showering and toileting. The chair can be wheeled over a standard toilet and used in a large shower alcove. Soft, adjustable, non-stick, quick-drying nylon back support Seat and back support easily detachable for cleaning or storage Swing away arm and foot supports Curved foot supports with side heel protection</p>	<p>Fully height adjustable from 450 - 600mm Seat height (at 20mm intervals) Maximum user weight 135kgs</p>	

No

Mobile Shower Commode,
Transit, Bariatric, 560mm wide,
Stainless Steel

Assists with the transportation of the user for showering and toileting. The chair can be wheeled over a standard toilet and used in a large shower alcove.

Wider seat widths and reinforced frame to suit the bariatric client

Swingback arms with padded armpads assist side transfers and can be positioned at two different heights

Adjustable safety arms for additional security

Sliding footrest is textured and weight-bearing for added safety

Footrest slides back under frame to assist patient transfers

PVC slatted backrest

Mounted on 4 x 125mm stainless steel locking castors

Durable stainless steel frame

Tested and conforms to Australian Standards AS/NZS 3973: 2009I

Seat Width - 560mm

Seat Depth - 430mm

Seat Height - 595mm


Armrest Height -
200/240mm



Backrest Height - 340mm

Overall Width - 640mm

Max User Weight - 300kg



No	Mobile Shower Commode, Transit, Bariatric 610mm wide, Stainless Steel	<p>Assists with the transportation of the user for showering and toileting. The chair can be wheeled over a standard toilet and used in a large shower alcove. Wider seat widths and reinforced frame to suit the bariatric client</p> <p>Swingback arms with padded armpads assist side transfers and can be positioned at two different heights</p> <p>Adjustable safety arms for additional security</p> <p>Sliding footrest is textured and weight-bearing for added safety</p> <p>Footrest slides back under frame to assist patient transfers</p> <p>PVC slatted backrest</p> <p>Mounted on 4 x 125mm stainless steel locking castors</p> <p>Durable stainless steel frame</p> <p>Tested and conforms to Australian Standard AS/NZS 2072-2000</p>	<p>Seat Width - 610mm</p> <p>Seat Depth - 430mm</p> <p>Seat Height - 595mm</p> <p>Armrest Height - 200/240mm</p> <p>Backrest Height - 340mm</p> <p>Overall Width - 690mm</p> <p>Max User Weight - 300kg</p>	
----	---	---	---	---

No	Mobile Shower Commode, Self Propelled, SS, 460mm	<p>Assists with the transportation of the user for showering and toileting. The chair can be wheeled over a standard toilet and used in a large shower alcove. Swingback arms with padded armpads for side transfers</p> <p>Adjustable safety arms for additional security</p> <p>Armrests can be positioned at two different heights</p> <p>Removable PVC slatted backrest</p> <p>Adjustable height swingaway legrests or slide out weight-bearing footrest</p> <p>Mounted on 2 x 125mm stainless steel castors and 2 x 22 inch rear wheels</p> <p>Durable stainless steel frame.</p>	<p>Seat Width - 460mm Seat Depth - 430mm Seat Height - 595mm Armrest Height - 200/240mm Backrest Height - 340mm Overall Width - 660mm Max User Weight - 200kg</p>	
No	Mobile Shower Commode, Self Propelled, Adjustable Height	<p>Easy to maneuver, for the independent user, the 24-inch wheels are positioned to make it easy to maneuver the chair without compromising user safety or comfort. It is easy to transfer to and from the chair thanks to the large, unrestricted seat in front of the wheels. The brake handles are also positioned below the seat so that they do not get in the way during transfer, while remaining easy to reach. Front seat recess; fully open at rear Long, removable, swing away arm supports Curved footplate with side protection for heels; slides in under seat Footplate features full height and tilt adjustment.</p>	<p>Width - 69 cm</p> <p>Width between arm supports - 43.5 cm</p> <p>Seat height - 55 cm</p> <p>Seat width - 48 cm</p> <p>Seat depth - 44 cm</p> <p>Space below seat - without pan holder: 49.5 cm, with pan holder: 48 cm Opening - depth x width, 36 cm x 21,5 cm</p> <p>Weight - 20 kg</p> <p>Maximum user weight - 130 kg</p>	

No

Denmark Mobile Shower
Commode

The Denmark mobile shower/commode chair is attendant propelled and has height adjustable and angle adjustable footrests. The slip resistant armrests can be removed sideways providing an advantage for side way transfers. The Chair also has 100 mm castors, an ergonomic cutout seat for easy hygiene access and a soft flexible plastic backrest. The weight capacity is 150kgs.

Weight (incl. Footrests):
12.6kg
Capacity: 150kg
Seat Depth: 45cm
Seat Height: 44.5 - 56.5cm
(adjustable)
Width Between Armrests:
47.5cm
Max Width: 58cm
Length (without the
footrests): 60cm



No

Aquatec Ocean VIP Ergo Tilt
Shower Commode

The Aquatec Ocean Ergo Family of shower chair commodes have been designed with user comfort, independence and dignity in mind. Based on the well established Aquatec Ocean range, the Ocean Ergo family comes packed with a range of new features, designed to make personal care activities easier for both the user and carer.

SWL - 150kg.
Seat Height 475 - 600mm
Seat width 480mm
Overall width 565mm
Overall Depth 1000mm
(including footrests)



No

ETAC Swift Mobil-2 Tilt in Space
Mobile Shower Commode with
Pan

The Swift Mobil Tilt-2, has a unique tilt feature that gives increased security for the user and an ergonomic working position for the carer. It's -5° to +35° tilt function can be operated from either side of the chair.

This allows the carer to reassure the user when tilting.

SWL - 160kg
Height 500 - 650mm
Width 500mm
Width 640mm

Seat
Seat
Overall



Shower Chairs

Yes

Shower Chair, Standard,
Adjustable - With Arms

A standard no-fuss, safe and reliable chair which is height adjustable. It features a contoured, plastic moulded seat, which is both comfortable & hygienic. It has been cleverly designed with a gentle slope, so that water can flow away for easy drainage. The steel frame using 22.2mm tubing enables users to gain a firm grip making transferring safer and easier. Tested to Australian standards.

Overall Height : 775mm - 900mm

Overall Width : 545mm

Overall Depth : 515mm

Seat Width: 440mm

Seat Depth: 405mm

Seat Height Range:

410mm–535mm

Width between Armrests:
490mm

Armrest (Seat to top of
handle): 200mm

Backrest Height: 350mm

Frame: Steel

Frame Diameter (armrest
and leg): 22.2mm

Extension Leg Diameter:
25.8mm

Extension Leg
(outer/inner): Outer

Rubber Tip (White): DIA
40mm x DIA 23.8mm



Colour: White



Net Weight: 3.7kg



Maximum User Weight:
125kg




Warranty: 12 months



Yes	Shower Chair - Bariatric - 525 Wide	<p>Reinforced wide frame – Provides stability for larger users</p> <p>Contoured plastic seat with drainage holes and windows - Provides postural support, improves drainage and reduces skin adhesion</p> <p>Height adjustable legs – Easy push button height adjustment for changing needs</p> <p>Additional height increments on one leg – Accommodates uneven surfaces</p> <p>All frames are individually treated for corrosion resistance</p> <p>Raised armrests – Assist the user transferring on and off the elevated seat</p> <p>Durable non-slip feet – Provide excellent grip and stability in wet conditions and a stainless steel washer within enhances durability</p>	<p>SWL 200kg</p> <p>Overall Height: 735 - 910 mm</p> <p>Overall Width: 570 mm</p> <p>Overall Depth: 570 mm</p> <p>Width Between Arms: 525 mm</p>	
Yes	Shower Chair - Bariatric - 550 Wide	<p>Height adjustable</p> <p>Wider seat widths and reinforced frame to suit the bariatric client</p> <p>Padded seat and curved backrest provide comfort and postural support</p> <p>Non-slip rubber feet</p> <p>Powder coated steel frame</p> <p>Electrophoretic surface treatment to prevent rusting</p>	<p>Seat Width 550mm</p> <p>Seat Height 470 - 640mm</p> <p>Base Width 650mm</p> <p>Base Depth 520mm</p> <p>Max User weight 300kg</p>	

Yes	Shower Chair - Bariatric - 600 Wide	<p>Height adjustable</p> <p>Wider seat widths and reinforced frame to suit the bariatric client</p> <p>Padded seat and curved backrest provide comfort and postural support</p> <p>Non-slip rubber feet</p> <p>Powder coated steel frame</p> <p>Electrophoretic surface treatment to prevent rusting</p>	<p>Seat Width 600mm</p> <p>Seat Height 470 - 640mm</p> <p>Base Width 700mm</p> <p>Base Depth 520mm</p> <p>Max User weight 300kg</p>	
Yes	Shower Chair - Bariatric - 610 Wide	<p>Backrest can be adjusted in depth +/- 60mm and has curved tubing to ensure comfort for any body shape.</p> <p>Having the right height for the user ensures success in getting in and out of the chair and these chairs are adjustable from 420 to 550mm.</p> <p>Butterfly armrests ensure the body doesn't get squeezed getting in to and out of the chair whilst providing a good grip for sitting and getting up.</p> <p>The high-tech, lightweight tubing provides the highest strength at the lowest weight</p> <p>Removable seat pad ensures these items are easy to clean.</p>	<p>Seat Width: 610mm</p> <p>Seat Depth: 400 to 520mm</p> <p>Overall Width: 730mm</p> <p>Weight: 7kg</p> <p>Maximum User Weight: 325kg</p>	

Yes	Shower Chair - Bariatric - 650 Wide	<p>Height adjustable</p> <p>Wider seat widths and reinforced frame to suit the bariatric client</p> <p>Padded seat and curved backrest provide comfort and postural support</p> <p>Non-slip rubber feet</p> <p>Powder coated steel frame</p> <p>Electrophoretic surface treatment to prevent rusting</p>	<p>Seat Width 650mm</p> <p>Seat Height 470-640mm</p> <p>Base Width 750mm</p> <p>Base Depth 520mm</p> <p>Max User Weight 300kg</p>	
Yes	Shower Chair, Commode & Toilet Aid, Bariatric (All in One)	<p>Redesigned with the carer in mind, the Bariatric All In One is an over toilet aid, shower chair and commode. Comes with square bowl with easy front access to remove bowl. The square bowl does not have a lid. Optional splash guard available. The seat lifts up and comes with a seat cover. The removable backrest adds stability when used as a bedside commode and can be removed when using over a toilet. The moulded armrests assist gripping for safe and easy transferring.</p>	<p>Total Width: 790mm</p> <p>Total Depth: 530mm</p> <p>Seat Depth: 450mm</p> <p>Seat Height: 465 - 615mm (to over seat top)</p> <p>Width Between Arms: 600mm</p> <p>Height Adjustable: 430-570mm</p> <p>Safe Working Limit: 295kg</p>	

Yes	Shower Chair, Swivel	<p>Height adjustable</p> <p>Swivel motion assists easy transfers from standing to showering position</p> <p>Seat locks in both transfer and showering positions</p> <p>Integrated arms assist transfers</p> <p>Drainage slots drain water away from seating area</p> <p>Non-slip rubber tips</p> <p>Powder coated steel frame</p>	<p>Width Between Arms 480mm</p> <p>Seat Height 480-560mm</p> <p>Base Width 560mm</p> <p>Base Depth 540mm</p> <p>Max User Weight 160kg</p>	
	Shower Stools			
	Shower Stool Adjustable - With Arms	<p>The standard shower stool is height adjustable and features side handles to assist with transfers. The stool has a small footprint and is ideal for smaller shower alcoves.</p>	<p>Overall Width : 510mm</p> <p>Overall Depth : 410mm</p> <p>Width Between Arms: 450mm</p> <p>Weight of Stool: 3.2kg</p> <p>Height Adjustable: 490 - 610mm</p> <p>Maximum User Weight: 110kg</p>	
Yes	Shower Stool Adjustable - Without Arms	<p>Height adjustable</p> <p>Side handles to assist transfers</p> <p>Durable plastic moulded seat</p> <p>Non-slip rubber tips</p> <p>Electrophoretic surface treatment to prevent rusting</p> <p>Powder coated steel frame</p>	<p>Width 380mm</p> <p>Seat height 510 - 660mm</p> <p>Base width 550mm</p> <p>Base Depth 380mm</p> <p>Max user weight 125kg</p>	

Yes

Shower Stool Adjustable - With Backrest

Height adjustable Padded seat and back for added comfort Integrated arms assist transfers Non-slip rubber tips Electrophoretic surface treatment to prevent rusting Powder coated steel frame

Width Between arms 460mm Seat Height 520-670mm Base Width 550mm Base Depth 380mm Max User Weight 125kg






Yes




Shower Stool Adjustable -
Heavy Duty Padded - Aluminium





Constructed from sturdy no rust
aluminium with easy push button
adjustment legs
The internal sliding legs provide
smooth and effortless height
adjustment whilst avoiding any chance
of water filling the legs
Reinforced frame using 25.4mm tubing
enables user to gain a firmer grip
making transferring safer & easier
The seat is padded for comfort and is
removable for cleaning

Overall Height (footprint):
679mm - 829mm
Overall Width (footprint):
655mm
Overall Depth (footprint):
445mm
Seat Width: 445mm
Seat Depth: 275mm
Seat Height Range:
539mm–689mm
Width between Armrests:
538mm
Armrest (Seat to top of
handle): 140mm
Frame: Aluminium
Frame Diameter (armrest
and leg): 25.4mm
Extension Leg Diameter:
28.6mm
Extension Leg
(outer/inner): Outer
Rubber Tip (White): DIA
40mm x DIA 27mm
Colour: White
Net Weight: 3.0kg
Maximum User Weight:
150kg
Warranty: 12 months
ARTG ID: 291268
Tested to Australian





Yes	Shower Stool With Removable Back and Arms - Bariatric	<p>Aluminum construction in light weight and corrosion free</p> <p>Cross-brace to enhance the stability during sitting</p> <p>Blow molded bench and back provide comfort and strength</p> <p>Drainage in seat to reduce slipping</p> <p>Limited Lifetime Warranty</p>	<p>Seat width 51cm</p> <p>Seat depth 36cm</p> <p>Seat height 36-51cm</p> <p>Overall width 58cm</p> <p>Overall depth 38cm</p> <p>Overall height 68-84cm</p> <p>Weight Capacity 227kg</p> <p>Total weight 3.3kg</p>	
	Accessories			
	Shower Hose Single Tap 2 Metre - Push on	<p>Fixes to existing tap spouts Lightweight plastic handpiece Supplied with two rubber reducing collars to suit most tap spouts Easily fitted and removed</p>	Length 2m	
Yes	Shower Hose Double Tap 2 Metre - Push on	<p>Hand Shower – Dual Tap Clamp-on. Fixes to existing tap spouts Lightweight plastic handpiece Supplied with two rubber reducing collars to suit most tap spouts Easily fitted and removed</p> <p>Length 2m</p>	Length 2m	
Toileting Aids				

Yes	Toilet Seat Extension - With Lid 50mm	- Ergonomic contoured one piece seat with lid Tapered flange with front and rear access Stain and odour resistant Easy to secure with two hand adjustable brackets	Height 50mm. Maximum user weight 130kg	
Yes	Toilet Seat Extension - With Lid 100mm	- Ergonomic contoured one piece seat with lid Tapered flange with front and rear access Stain and odour resistant Easy to secure with two hand adjustable brackets	Height 100mm. Maximum user weight 130kg	
Yes	Toilet Seat Extension - 50mm	Increases toilet seat height Ideal for people with hip or knee injuries or those who require greater seat height Easily positioned and removed using location lip and adjustable clamps Front and rear cutouts for personal cleaning Contoured seat for greater comfort Moulded plastic is stain and odour resistant	Height: 50mm Max user weight: 190kg	

Yes	Toilet Seat Extension - 100mm	Increases toilet seat height Ideal for people with hip or knee injuries or those who require greater seat height Easily positioned and removed using location lip and adjustable clamps Front and rear cutouts for personal cleaning Contoured seat for greater comfort Moulded plastic is stain and odour resistant	Height: 100mm Max user weight: 190kg	
Yes	Toilet Surround - Adjustable	Frame to assist toilet transfers Height adjustable Non-slip rubber tips Powder coated steel	Height - 625-775 Width between Arms - 455 Base Width - 540 Base Depth - 490 Max User Weight - 125kg	
Yes	Toilet Surround - Wide	Heavy duty fixed height frame to assist toilet transfers Fitted with rear cross bar for additional strength and stability Non-slip rubber tips Powder coated steel	Height - 685 Width between Arms - 620 Base Width - 670 Base Depth - 520 Max User Weight - 130kg	
Yes	Toilet Seat Raiser - Adjustable with Splash Guard (SWL 125kg)	The Over Toilet Aid assists transfers on & off the toilet. This toilet surround frame is height adjustable, with a push button adjustment. Powder coated steel frame with non slip rubber feet. Comes complete with a seat and splashguard	Seat Height 440-590mm Seat Width 370mm Width Between Arms 450mm Max User Weight 125kg	

Yes	Toilet Seat Raiser 525 - BARIATRIC Adjustable with Splash Guard (SWL 175kg)	<p>Raised armrests – Assist the user transferring on and off the Over Toilet Aid</p> <p>Additional height increments on one leg – Accommodates uneven surfaces</p> <p>Height adjustable legs – Easy push button height adjustment for changing needs</p>	<p>Seat Height 450 - 630mm</p> <p>Seat Depth 435mm</p> <p>Width Between Arms 525mm</p> <p>Max User Weight 175kg</p>
Yes	Toilet Seat Raiser - Bariatric 550mm Wide	<p>Provides a raised seat to assist in toileting and transfers.</p> <p>Wider seat widths and reinforced frame to suit the bariatric client.</p> <p>Ideal for people with hip or knee injuries or those who require greater seat height</p> <p>Ergonomic angled arms assist transfers</p> <p>Supplied with lift-back toilet seat and removable splashguard</p> <p>Non-slip rubber tips</p> <p>Powder coated steel</p>	<p>Width Between Arms 550mm</p> <p>Seat height 550mm</p> <p>Armrest height 230mm</p> <p>Base width 680mm</p> <p>Base Depth 610mm</p> <p>Maximum user weight 200kg</p>



Yes	Toilet Seat Raiser - Bariatric 600mm Wide	<p>Provides a raised seat to assist in toileting and transfers</p> <p>Wider seat widths and reinforced frame to suit the bariatric client</p> <p>Ideal for people with hip or knee injuries or those who require greater seat height</p> <p>Ergonomic angled arms assist transfers</p> <p>Supplied with lift-back toilet seat and removable splashguard</p> <p>Non-slip rubber tips</p> <p>Powder coated steel</p>	<p>Width Between Arms 600mm</p> <p>Seat height 550mm</p> <p>Armrest height 230mm</p> <p>Base width 730mm</p> <p>Base Depth 610mm</p> <p>Maximum user weight 250kg</p>	
Yes	Toilet Seat Raiser - Bariatric 650mm Wide	<p>Provides a raised seat to assist in toileting and transfers</p> <p>Wider seat widths and reinforced frame to suit the bariatric client</p> <p>Ideal for people with hip or knee injuries or those who require greater seat height</p> <p>Ergonomic angled arms assist transfers</p> <p>Supplied with lift-back toilet seat and removable splashguard</p> <p>Non-slip rubber tips</p> <p>Powder coated steel</p>	<p>Width Between Arms 650mm</p> <p>Seat height 550mm</p> <p>Armrest height 230mm</p> <p>Base width 780mm</p> <p>Base Depth 610mm</p> <p>Maximum user weight 300kg</p>	

Yes

Shower Chair, Commode &
Toilet Aid, Bariatric (All in One)

The bariatric All in One provides for being a toilet seat raise, shower chair and commode chair. The legs are height adjustable. It has a square bowl with easy access for removal from the front of the chair. The square bowl does not have a lid. The seat lifts up and comes with a seat cover. The removable backrest adds stability when used as a bedside commode and can be removed when using over a toilet. The moulded armrests assist gripping for safe and easy transferring.

Total Width: 790mm
Total Depth: 530mm
Seat Depth: 450mm
Seat Height: 465 - 615mm
(to over seat top)
Width Between Arms:
600mm
Height Adjustable: 430-
570mm
Safe Working Limit: 295kg



CHSP EQUIP IN
SCOPE LIST

PRODUCT

PRODUCT DESCRIPTION

SPECIFICATIONS

IMAGE

BEDROOM

Bedside Commode

Yes

Commode - Adjustable Height

The Bedside Commode has a clip on toilet seat, which clips on in seconds. There are no metal fittings on toilet seat. It is also height adjustable and easy to clean. It has a light blue deluxe padded overseat, as well as a removable round bowl with lid and handle. The lid eliminate spillages and reduces risks of cross infection. It also features a powder coated dark blue frame. This bedside commode looks very smart in a bedroom.

Overall Height: 820 - 970mm
Overall Width: 610mm
Overall Length: 550mm
Seat Height: 460 - 610mm
Weight: 8kg
Maximum User Weight: 110kg





Yes

Classic Bedside Commode

Allows for toileting in the bed room for users unable to travel between the bedroom and bathroom
Available in Champagne or Slate Vinyl
Raised Padded Armrests
Easy Clean Vinyl Surfaces
Height Adjustable
Removable, Padded Seat
Non-Slip Rubber Tips

Seat Height: 475 - 615 mm



Yes	Commode - Adjustable Height - Heavy Duty - 580mm	Steel Heavy Duty Commode With removable Backrest. Durable and rigid frame in extra wide for high weight loading Adjustable height.	Seat 58cm width x 46cm depth Max User Weight 200kg.	
	Commode - Bariatric - 650mm Wide	Wider seat widths and reinforced frame to suit the bariatric client Padded seat, padded armrests, lift-back toilet seat and removable commode bowl Adjustable seat height Padded seat cover slides back easily for access to commode Three sided vanity skirt conceals commode bowl Powder coated steel frame	Width between arms 650mm Seat height: 390 - 540mm Armrest Height: 300mm Max user weight 150kg	

Yes

Discuss with CL

Shower Chair, Commode & Toilet Aid, Barriatric (All in One)

Redesigned with the carer in mind, the Bariatric All In One is an over toilet aid, shower chair and commode. Comes with square bowl with easy front access to remove bowl. The square bowl does not have a lid. Optional splash guard available. The seat lifts up and comes with a seat cover. The removable backrest adds stability when used as a bedside commode and can be removed when using over a toilet. The moulded armrests assist gripping for safe and easy transferring.

Total Width: 790mm
Total Depth: 530mm
Seat Depth: 450mm
Seat Height: 465 - 615mm (to over seat top)
Width Between Arms: 600mm
Height Adjustable: 430-570mm
Safe Working Limit: 295kg



Bed Sticks



Horizontal Bed Rail




This Bed Rail is designed to prevent the user from rolling out of bed. Once the user is in bed, simply slide the Bed Rail under the mattress. The weight of the user secures the Bed Rail. Constructed from sturdy epoxy coated steel.




Height: 500mm
Width (Base): 580mm
Length (Rail): 600mm
Weight: 2.5kg



Bed Positioning

Yes	Bed Bolster - Mattress Extension	<p>200 mm long foam mattress extension for additional pressure relief under heel zone, available in Standard and King Single widths.</p> <p>Suitable for use in most standard residential care, hospital and home care beds with sleeping platform extension capabilities. Standard model suits any 150 mm high foam mattress; King Single model suits Triple Layer King surface.</p>	<p>King Single (180 high, 1050mm wide)</p> <p>Single (150 high, 865mm wide)</p> <p>Features</p> <ul style="list-style-type: none">high-density memory foam upper layer with medium density foam core basedouble-strength multi-stretch waterproof cover reduces the potential for shear injuriesseams are stitched, welded and taped for maximum infection control <p>Standard and King Single widths available</p>	
	Back Rest - Adjustable	<p>Multi-angled pillow support</p> <p>Provides assistance and support in bed</p> <p>Wide aluminium frame fitted with strong breathable nylon back</p> <p>Angle can be adjusted via a simple ratchet mechanism</p>	<p>Length 495mm</p> <p>Width 700mm</p> <p>Height 490mm</p>	

Yes	Foam Bed Wedge	<p>Angled foam cushion provides comfort and support to the upper body</p> <p>Ideal for those who suffer from respiratory and reflux problems</p> <p>Removable quilted cotton cover</p> <p>Wedge can be turned to provide two different recline angles (15 degrees or 75 degrees)</p> <p>Foam insert can be removed to alter the angle by approx 5 degrees</p>	<p>Length 680mm</p> <p>Width 710mm</p> <p>Height 240-600mm</p>	
	Bed / Chair Blocks			
	BED Blocks Plastic 40mm	<p>Plastic block raisers designed to increase the height of a bed.</p> <p>Increases height by 40 mm</p> <p>Stackable</p>	<p>Internal radius: 65 mm</p> <p>Suits square legs up to 45 mm x 45 mm</p> <p>Weight Limit: Furniture plus user must weigh less than 250 kg</p>	
Yes	CHAIR Blocks Plastic 40mm	<p>Plastic block raisers designed to increase the height of a chair.</p> <p>Increases height by 40 mm</p> <p>Stackable</p>	<p>Internal radius: 65 mm</p> <p>Suits square legs up to 45 mm x 45 mm</p> <p>Weight Limit: Furniture plus user must weigh less than 250 kg</p>	

Yes	BED Blocks Plastic 100mm	Plastic block raisers designed to increase the height of a bed. Increases height by 100 mm Stackable	Internal radius: 65 mm Suits square legs up to 45 mm x 45 mm Weight Limit: Furniture plus user must weigh less than 250 kg	
	CHAIR Blocks Plastic 100mm	Plastic block raisers designed to increase the height of a chair. Increases height by 100 mm Stackable	Internal radius: 65 mm Suits square legs up to 45 mm x 45 mm Weight Limit: Furniture plus user must weigh less than 250 kg	
	Overbed Tables			
Yes	Overbed Table	Eoxy coated frame with a vinyl top. Recessed edges, which enables the table to be cleaned easily. The 50mm twin wheel castors allow the table to be easily moved away from the bed when not in use.	Length 800mm Overall Width 400mm Minimum Height 710mm Maximum Height 1050mm Height of Castor Wheels 50mm	

Yes

Overbed Table - Split Top Tilt

Recessed Thermoform top to prevent spills

Lightweight mobile tilting table positions over bed or chair

Can be used for eating, drawing or other activities

Angle adjustable to 45 degrees, ideal for homecare environment

Side table does not tilt, suitable for drink and food

Fully height adjustable

"C" base shape fitted with four castors enables close positioning to the bedside or chair

Stylish cream colour finish

Powder coated mild steel frame construction

Height: 780 - 1140mm

Width: 430mm

Length: 800mm



Self Help Poles/Ladders

Yes

Bed Cradle - Adjustable Height

The bed cradle is designed to keep bed linen off sensitive areas of the patient's lower leg or feet by creating clearance between bed linen and the mattress. The lower frame slides underneath the mattress, securing the cradle to the bed, while the upper frame is positioned horizontally over the lower limbs. Bedding is then placed over the frame to remove any weight or pressure from the sensitive limbs below. Height adjustable frame allows the user to customise the space between mattress surface and bedding.

Height adjustable 440 - 640 mm



Yes

Bed Ladder

A solution for clients who have physical difficulties sitting up in or getting out of bed. Step by step they gradually pull themselves up from a lying to a sitting position. The action involved also strengthens the stomach and arm muscles of the user. Can also be used to assist in transfers from wheelchairs to beds. Plastic ladder rungs give a firm and safe grip. Cord ends tie around the bed legs. Suitable for single or double beds.

Rung Width 195mm
Rung Diameter 25mm
Cord Length 3m (to first rung)



Yes

Self Help Pole - Freestanding

Free standing support handle provides assistance for the user to reposition or transfer in and out of bed

Handle height is easily adjusted to suit the user

Sits behind or beside the bed

Safe work load 150kg



CHSP EQUIP
IN SCOPE LIST

PRODUCT

PRODUCT DESCRIPTION

SPECIFICATIONS

IMAGE

MOBILITY AIDS

Rollator Frames

Yes

Rollator Walking Frame - Mini

Provides support to assist a person to walk who has difficulty due to balance, pain or injury. Frames comes set up with 2 wheels at the front and 2 skis at the back for easier maneuvering. Frame folds up flat for transporting.

Max height 82cm
Min height 72cm.
Max user weight. 150kg



Yes

Rollator Walking Frame - Standard 5 inch wheels

Provides support to assist a person to walk who has difficulty due to balance, pain or injury. Frames comes set up with 2 wheels at the front and 2 skis at the back for easier maneuvering. Frame folds up flat for transporting.

Adjustable height 810 - 910mm.
Overall width 670mm.
Width between handles 430mm
Depth 500mm.
Recommended mainly for indoor use
Max user weight. 150kg



Yes

Rollator Walking Frame - Bariatric 227kgs

Oversize frame, wider and deeper to accommodate individuals up to 227kgs
Provides support to assist a person to walk who has difficulty due to balance, pain or injury. Frames comes set up with 2 wheels at the front and 2 skis at the back for easier maneuvering. Frame folds up flat for transporting.

Maximum user weight: 227kg
Height adjustable: 32.5"~39.5" with 5 wheels



Pick Up Walking Frames

Yes

Pick Up Walking Frame - Mini

Height adjustable
Folds flat for easy storage and transportation
Foam handles
Lightweight aluminium
Tip size 28mm

Handle Height 700 - 805mm
Max Base Size 630 x 500mm
Weight 2.5kg
Max user weight 150kg



Yes

Pick Up Walking Frame -- Standard

Folding aluminium frame with foam rubber hand grips.

Adjustable height 770 - 940mm with rubber ferrules on all 4 legs.
Width 455mm, Depth 400mm.



Yes



Pick-up Walking Frame - Bariatric 227kgs



Oversize frame design provides wider and deeper to accommodate individuals up to 227kgs
Precision design provides additional strength while adding minimal weight



Maximum user weight 227kgs
Height adjustable: 32.5"~39.5"








Walking Sticks

Yes	Walking Sticks - T Handle	<p>T shaped handle, single point, adjustable. To provide support people who have difficulty walking due to balance problems, pain or injury.</p> <p>The height of the stick can be adjusted to match the height of the user an important feature when selecting the stick that is right for you. Adjust the height of the telescopic leg with the push of a button. The tip of the stick is covered in a solid rubber tip, which can be replaced should it show signs of wear after prolonged use.</p>	<p>Handle height: 79 cm - 102 cm (25 mm increments) Approximate user height: 5' 3" - 6' 1" Rubber tip/ferrule: 19 mm Product weight: 300 gm Maximum user weight capacity: 114 kg</p>	
Yes	Walking Stick - Bariatric - 225 Kg	<p>This sturdy steel walking stick provides support to a person who has difficulty walking due to balance problems, pain or injury.</p> <p>The height of the stick can be adjusted to match the height of the user an important feature when selecting the stick that is right for you. Adjust the height of the telescopic leg with the push of a button. The tip of the stick is covered in a solid rubber tip, which can be replaced should it show signs of wear after prolonged use.</p>	<p>Rubber Hand Grip Handle Height: 73 cm - 96 cm (25 mm Increments) Approximate User Height: 5' 2" - 6' 1" Rubber tip/ferrule: 22 mm Maximum user weight capacity: 225 kg</p>	

Yes	Walking Stick - Crook Handle	<p>To provide support people who have difficulty walking due to balance problems, pain or injury. Constructed from aluminium to allow durability and strength whilst being light in weight.</p> <p>The height of the stick can be adjusted to match the height of the user an important feature when selecting the stick that is right for you. Adjust the height of the telescopic leg with the push of a button.</p> <p>The tip of the stick is covered in a solid rubber tip, which can be replaced should it show signs of wear after prolonged use.</p>	<p>Height adjustment range: 760mm - 990mm</p>	
Yes	Walking Stick - Ergonomic - Left Hand	<p>Anatomical Walking StickConstructed from aluminium to allow durability and strength whilst being light in weight.</p> <p>Supplied with anti-slip rubber tip. This walking stick is available in both a left hand and right hand option to ergonomically fit the natural shape of the user's hand for maximum comfort and ease of use.</p>	<p>Height adjustment range: 760mm - 990mm</p>	

Yes	Walking Stick - Ergonomic - Right Hand	<p>Anatomical Walking StickConstructed from aluminium to allow durability and strength whilst being light in weight. Supplied with anti-slip rubber tip. This walking stick is available in both a left hand and right hand option to ergonomically fit the natural shape of the user's hand for maximum comfort and ease of use.</p>	:Height adjustment range: 760mm - 990mm	
Yes	Walking Stick - Folding	<p>A lightweight, aluminium walking stick with compact folding mechanism for convenient storage. The stick folds and unfolds simply and easily, making it ideal for everyday use or travel.</p> <p>The height of the stick can be adjusted to match the height of the user an important feature when selecting the stick that is right for you. Adjust the height of the telescopic leg with the push of a button. The tip of the stick is covered in a solid rubber tip, which can be replaced should it show signs of wear after prolonged use.</p>	<p>lightweight folding walking stick height adjustable easy to collapse and unfold maximum user weight 113kg black colour</p>	

Yes	Walking Stick - (Swan neck) with strap	<p>To provide support people who have difficulty walking due to balance problems, pain or injury. Constructed from aluminium to allow durability and strength whilst being light in weight.</p> <p>The height of the stick can be adjusted to match the height of the user an important feature when selecting the stick that is right for you. Adjust the height of the telescopic leg with the push of a button. The tip of the stick is covered in a solid rubber tip, which can be replaced should it show signs of wear after prolonged use.</p>	<p>Height adjustment range: 810mm - 1035mm</p> <p>Tube size: Upper tube 22.2mm, lower tube (telescopic) 19mm</p> <p>Rubber tip (Outer size Ø38XH48; Inner sizeØ18.5xD34)</p> <p>Maximum User Weight: 113kg</p> <p>Colour: Silver</p>	
Yes	Walking Stick - (Swan neck) without strap	<p>To provide support people who have difficulty walking due to balance problems, pain or injury. Constructed from aluminium to allow durability and strength whilst being light in weight.</p> <p>The height of the stick can be adjusted to match the height of the user an important feature when selecting the stick that is right for you. Adjust the height of the telescopic leg with the push of a button. The tip of the stick is covered in a solid rubber tip, which can be replaced should it show signs of wear after prolonged use.</p>	<p>Height adjustment range: 810mm - 1035mm</p> <p>Tube size: Upper tube 22.2mm, lower tube (telescopic) 19mm</p> <p>Rubber tip (Outer size Ø38XH48; Inner sizeØ18.5xD34)</p> <p>Maximum User Weight: 113kg</p> <p>Colour: Silver</p>	
Quad Sticks				

Yes	Quad Stick Small - 15cm x 20cm	This lightweight aluminium quad base walking stick is suitable for a user with limited mobility who requires a four point base for extra stability. The offset foam handle can rotate for left or right hand use. It is adjustable in height, wrist strap provided,	Handle height: 74 cm - 97 cm (25 mm increments) Soft foam handgrip Base: 20 cm x 15 cm Rubber tips/ferrule: 13 mm Product weight: 800 gm Maximum user weight: 114 kg	
	Quad Stick Large - 17cm x 28cm	This lightweight aluminium quad base walking stick is suitable for a user with limited mobility who requires a four point base for extra stability. The offset foam handle can rotate for left or right hand use. It is adjustable in height, wrist strap provided,	Handle height: 76 cm - 98 cm (25 mm increments) Soft foam handgrip Base: 28 cm x 17 cm Rubber tips/ferrule: 16 mm Product weight: 800 gm Maximum user weight capacity: 113 kg	
	Quad Stick - Off Set - 23 x 39cm	This lightweight aluminium quad base walking stick is suitable for a user with limited mobility who requires a four point base for extra stability. The offset foam handle can rotate for left or right hand use. It is adjustable in height, wrist strap provided, Adjustable height, larger base.	Height 795-1020mm Base width 230mm Base Depth 390mm Tip size 16mm SWL 75kg	
4 Wheel Walkers				

Yes

4 Wheel Walker - CAREQUIP
MINI

Frame and backrest fold compact for storage and transportation
Height adjustable handles
Narrow frame and low seat height ideal for the smaller client, or where space is limited
Handbrakes can be locked for added stability
Supplied with removable backrest, vinyl bag and padded seat
Walker can be folded with vinyl bag still attached
Curved padded backrest for additional comfort and support
*Note: Product comes with a vinyl bag and not a basket as pictured.

Seat height 470mm Seat width 310mm Max user weight 130kg
Unit weight 6.2kg





Yes

4 Wheel Walker - STANDARD
BLUE

Frame and backrest fold compact for storage and transportation
Height adjustable handles
Soft-lock handbrakes
Curved padded backrest for additional comfort and support
Supplied with vinyl bag and padded seat
Walker can be folded with vinyl bag still attached

Seat height 525mm Seat width 350mm Max user weight 130kg
Unit weight 5.9kg



Yes	4 Wheel Walker - STANDARD RED	<p>Frame and backrest fold compact for storage and transportation</p> <p>Height adjustable handles</p> <p>Soft-lock handbrakes</p> <p>Curved padded backrest for additional comfort and support</p> <p>Supplied with vinyl bag and padded seat</p> <p>Walker can be folded with vinyl bag still attached</p>	<p>Seat height 525mm</p> <p>Seat width 355mm</p> <p>Max user weight 130kg</p> <p>Unit weight 5.9kg</p>	
Yes	4 Wheel Walker - DAYS RED	<p>The Days 100 Series aluminium frame ensures a combination of lightweight design and robust build quality with a maximum user weight of 125kg.</p> <p>6" wheels</p> <p>Ergonomic handgrips to make them easy to grip.</p> <p>Adjustable height handles with 5 settings at 25mm intervals</p> <p>Rest Seat</p> <p>Storage Bag</p> <p>Loop lockable brakes provide added security and convenience. The brake tension can also be easily adjusted to suit requirements</p> <p>Foldable backrest</p> <p>Folds easily for storage & transport - simply fold together by pulling the strap on the seat</p>	<p>Max User Weight: 125kg</p> <p>Seat Height: 500mm</p> <p>Seat Width: 350mm</p> <p>Seat Depth: 320mm</p> <p>Width Between Arms: 460mm</p> <p>Backrest Height: 750mm</p> <p>Overall Height: 755-875mm</p> <p>Overall Width: 620mm</p> <p>Overall Depth: 640mm</p> <p>Backrest: Straight, folding</p> <p>Wheel: 6"</p> <p>Tyre: PVC</p> <p>Net Weight: 6kg</p>	

Yes

DAYS Walker Boxed -
ADJUSTABLE Standard Red

Easy to assemble - no tools required
Retail carrying box with handle
Height adjustable seat and handles
Padded seat with storage compartment
bag underneath
Padded backrest
Hand operated brake system

Max User Weight: 136kg
with foldable backrest
Seat Height: 500 -
600mm
Seat Width: 356mm
Seat Depth: 305mm
Overall Height: 785 -
1000mm
Overall Width: 622mm
Overall Depth: 680 -
800mm
Width Between Handles:
460mm
Backrest: Curved, folding
Wheel: 7.5"
Tyre: PVC, Black
Net Weight: 7.4kg



Yes

DAYS Walker Boxed-
ADJUSTABLE Standard Blue

Easy to assemble - no tools required
Retail carrying box with handle
Height adjustable seat and handles
Padded seat with storage compartment
bag underneath
Padded backrest
Hand operated brake system

Max User Weight: 136kg
with foldable backrest
Seat Height: 500 -
600mm
Seat Width: 356mm
Seat Depth: 305mm
Overall Height: 785 -
1000mm
Overall Width: 622mm
Overall Depth: 680 -
800mm
Width Between Handles:
460mm
Backrest: Curved, folding
Wheel: 7.5"
Tyre: PVC, Black
Net Weight: 7.4kg



Yes

4 Wheel Walker - BLUE (LARGE WHEELS)

8" wheels
Ergonomic handgrips and adjustable height handles
Rest Seat
Storage Bag
Loop lockable brakes
Foldable backrest
Folds easily for storage & transport

Max User Weight: 125kg
Seat Height: 565mm
Seat Width: 350mm
Seat Depth: 320mm
Width Between Arms: 460mm
Backrest Height: 800mm
Overall Height: 840-940mm
Overall Width: 620mm
Overall Depth: 645mm
Backrest: Straight, folding
Wheel: 8"
Tyre: PVC
Net Weight: 7.2kg



Yes

4 Wheel Walker - GOLD (LARGE WHEELS)

8" wheels
Ergonomic handgrips and adjustable height handles
Rest Seat
Storage Bag
Loop lockable brakes
Foldable backrest
Folds easily for storage & transport

Max User Weight: 125kg
Seat Height: 565mm
Seat Width: 350mm
Seat Depth: 320mm
Width Between Arms: 460mm
Backrest Height: 800mm
Overall Height: 840-940mm
Overall Width: 620mm
Overall Depth: 645mm
Backrest: Straight, folding
Wheel: 8"
Tyre: PVC
Net Weight: 7.2kg



Yes

4 Wheel Walker - HEAVY DUTY

To provide support and stability whilst walking. Seat lifts up to allow for storage in basket below
7.5" wheels
Adjustable height handles
Rest seat
Storage bag
Loop lockable brakes
Folds easily for storage & transport
Screen reader support enabled.

Max User Weight: 200kg
Seat Height: 570mm
Seat Width: 495mm
Seat Depth: 355mm
Overall Height: 855 - 975mm
Overall Width: 735mm
Overall Depth: 655mm
Backrest: Curved
Wheel: 7.5"
Tyre: PVC, Black
Net Weight: 8.8kg
Screen reader support enabled.



Yes

4 Wheel Walker -PUSH DOWN BRAKES

To provide support and stability whilst walking. Seat lifts up to allow for storage in basket below. The unique brake system is activated by pushing down on the handles or whenever you sit on the seat. The robust yet lightweight alloy frame has tube-in-tube construction for increased strength and longer lifespan. The castors are made from solid rubber, making them free from punctures. They are designed for easy operation and good manoeuvrability especially on indoor surfaces. Their slightly smaller size makes them especially good for negotiating doorways or tight spaces. The orthopedic moulded handgrips are designed to contour entirely to your hand, minimising stress or strain to the wrist. The height of the handles is easy to adjust, with an extension range of 13 cm to ensure that the handles can be set to best match the height of the user.

With a compact fold down design, this walker is easy to transport or store when not in use.

Suitable for people up to 130 kilograms in weight.
Height 945mm max - 810mm min
Length 610mm
Width 575mm
Seat Depth 310mm
Seat Width 360mm
Seat height 557mm
Folded width 215mm
Castor size 153mm



Yes

Gutter Frame

Sturdy wheeled walking frame for indoor mobility training and rehabilitation. Designed for users who require the ability to weight bear through the forearms. Available in both non folding and folding models.

Note: Not designed to help users stand from a seated position

Forearm cuff height: 93 cm - 123 cm (25 mm increments)
Width: Front 40 cm / Rear 60 cm
Front swivel wheels and rear ski's for maneuverability
Adjustable hand grips to suit users forearm length
Maximum user weight capacity: 110 kg



Tri Walkers

Yes

Tri Walker with Tray and Basket

The Aluminium Tri Wheel Walkers have a lightweight aluminium frame that makes these walkers extremely easy to manoeuvre making them ideal for use indoors or out.

Ergonomic handgrips for optimum comfort
Adjustable height handles with 5 settings at 25mm intervals
Loop lockable brakes provide added security and convenience. The brake tension can also be easily adjusted to suit requirements
Vinyl bag supplied as standard
Folds easily for storage or transport using the simple easy to grab pull handle
Supplied with basket and tray.



Crutches




Yes



Crutches Axillary Large

Aluminium, push button adjustable height axilla crutches with rubber ferrules.

Maximum User Weight - 136kg
135 - 155cm



Yes	Crutches Axillary Medium	Aluminium, push button adjustable height axilla crutches with rubber ferrules.	Maximum User Weight - 136kg 114 - 135cm	
	Crutches Axillary Youth/Small	Durable and strong, these Aluminium Underarm Crutches are constructed from high tensile extruded aluminium making them both lightweight and very strong. Fitted with high quality rubber tips for stability. Colour Grey. Maximum User Weight - 136kg	Maximum User Weight - 136kg 94 - 114cm	
	Elbow Crutches Medium	Durable and strong, these Forearm Crutches are constructed from high tensile extruded aluminium, making them both lightweight and very strong. Fitted with high quality rubber tips for stability.	Maximum User Weight - 160kg 89 - 120cm	

Yes	Elbow Crutches Medium Ergonomic Grip	<p>Durable and strong, these Ergonomic Grip Forearm Crutches are constructed from high tensile extruded aluminium making them both lightweight and very strong. Ergonomic palm grip offers pressure distribution for maximum comfort. Fitted with high quality rubber tips for stability. Both hand to floor and palm grip to arm cuff sections are adjustable.</p>	<p>Maximum User Weight - 160kg Hand Grip Height - 66 - 94cm, Arm Ring Adjustment- 20.5 - 27.5cm</p>	
Yes	Elbow Crutches Large Ergonomic Grip	<p>Durable and strong, these Ergonomic Grip Forearm Crutches are constructed from high tensile extruded aluminium making them both lightweight and very strong. Ergonomic palm grip offers pressure distribution for maximum comfort. Fitted with high quality rubber tips for stability. Both hand to floor and palm grip to arm cuff sections are adjustable.</p>	<p>Maximum User Weight - 160kg Hand Grip Height - 83 - 111cm, Arm Ring Adjustment- 20.5 - 27.5cm</p>	
Electric Scooters				

No

745 Plus

Stylish eye catching design with attractive modern colour range

This new design offers generous foot space & leg room, also with large battery capacity

Loop handlebars offer comfortable finger or thumb throttle operation

Adjustable high back contoured seat with simple operation of slide and swivel mechanisms to allow maximum user comfort

Infinite adjustment of steering tiller allowing everyone to easily find a comfortable position

High mount charger socket means easy access when charging

11" wheels which provide extra ground clearance

Adjustable suspension to give maximum comfort

4 pole motor delivering more power

Adjustable power saving LED headlight, indicators front and rear, rear tail lights and brake lights

Built in speed reduction switch which is activated when turning offering exceptional safety.

Max speed: 10kph

Range: 30-40km

Weight: 98kg

SWL: 150kg



No

Fende

Stylish new modern design with a modesty panel, giving this scooter a smart and attractive look

Adjustable prestige car-like front & rear suspension delivers a smooth ride over varied terrain

Comfortable and contoured leather feel high back seat with slide, swivel, recline and adjustable armrests

Unique battery storage design, allowing for greater floor space and a secure lockable storage area for personal items

New modern LCD dashboard display, making reading the battery, speed, time & temperature gauges clear and easy

LED lighting system delivers brighter illumination with less power consumption

Built-in speed reduction switch, which is activated when turning offering exceptional safety

With new easy to access 12 volt socket (for phone or GPS)

Max speed: 10kph

Range: 30km

Weight: 114.5kg

SWL: 136kg



No

Pegasus Pro

Comfort, safety and reliability with added performance. Ideal for longer trips into town or afternoons driving around the park with your grandchildren, the new ergonomically designed Pegasus PRO is a powerful upgrade to our standard Pegasus METRO model.

Max speed: 10kph
Range: 52km
Weight: 136kg
SWL: 160kg



CHSP EQUIP
IN SCOPE LIST

PRODUCT

PRODUCT DESCRIPTION

SPECIFICATIONS

IMAGE

WHEELCHAIRS

Self Propelled Wheelchair

Yes

Self-Propelled 14" SWIFT

The Days Swift Self-Propelled 14" Wheelchair has been designed to be a light, practical and easy to transport wheelchair, making it an ideal occasional use chair for both indoor and outdoor use. Its compact dimensions allow good manoeuvrability indoors & the light construction and ease to fold for storage or transport make it a very versatile chair.

Light frame 7.8kg
Folding backrest (for storage or transport)
Quick release wheels (for transport)
Fixed arms
Swing away and removable leg rests
Solid tyres (no maintenance required)
Extended brake handle (ease of brake application)
Maximum User Weight: 80kg



Yes

Self-Propelled 16" LINK
SPECIAL ORDER

Days Premium Link Wheelchair is a lightweight adaptable wheelchair suitable for full-time users and is designed for maximum practicality and adjustability.

Light weight aluminium - 17.2kg

- Quick release wheels
- Solid tyres
- 16" seat width
- Swing in/out legrests
- Swing back removable armrests
- Height adjustable back upholstery
- Adjustable castors
- Multi position rear wheels
- Anti-tipper comes as standard inclusion



Yes

Self-Propelled 16" SWIFT

The Days Swift Self-Propelled Wheelchair with Handbrakes has been designed to be a light, practical and easy to transport wheelchair, making it an ideal occasional use chair for both indoor and outdoor use. Its compact dimensions allow good manoeuvrability indoors & the light construction and ease to fold for storage or transport make it a very versatile chair.

Light frame 7.8kg

Folding backrest (for storage or transport)

Quick release wheels (for transport)

Fixed arms

Swing away and removable leg rests

Solid tyres (no maintenance required)

Handbrakes for added safety*

Extended brake handle (ease of brake application)

Maximum User Weight: 115kg



Yes

Self-Propelled 18" LINK

Days Premium Link Wheelchair is a lightweight adaptable wheelchair suitable for full-time users and is designed for maximum practicality and adjustability.

Light weight aluminium - 17.2kg

- Quick release wheels
 - Solid tyres
 - 18" seat width
 - Swing in/out legrests
 - Swing back removable armrests
 - Height adjustable back upholstery
 - Adjustable castors
 - Multi position rear wheels
- Can fit stump supports
- Drum brakes on the 20"-22"
 - Anti-tipper comes as standard inclusion



Yes

Self-Propelled 18" SWIFT

The Days Swift Self-Propelled Wheelchair with Handbrakes has been designed to be a light, practical and easy to transport wheelchair, making it an ideal occasional use chair for both indoor and outdoor use. Its compact dimensions allow good manoeuvrability indoors & the light construction and ease to fold for storage or transport make it a very versatile chair.

Light frame 7.8kg
Folding backrest (for storage or transport)
Quick release wheels (for transport)
Fixed arms
Swing away and removable leg rests
Solid tyres (no maintenance required)
Handbrakes for added safety*
Extended brake handle (ease of brake application)
Maximum User Weight: 115kg



Yes

Self-Propelled 18" WHIRL

The Healthcare Whirl is a versatile wheelchair with a range of features to accommodate different user needs. This sleek and stylish wheelchair is practical and will look good when you are out and about. Features: Quick release wheels (for transport) Removable armrests Swing away and removable leg rests Solid tyres (no maintenance required)

Specifications: * Weight: 16.8kg (complete) * User Weight: 120kg * Seat Width: 450mm * Seat Depth: 430mm * Backrest Height: 440mm * Seat Height - Front: 500mm * Seat Height - Rear: 475mm * Front Castor Wheel: 8" x 1" PU * Rear Wheel: 24" x 1 3/8" PU * Overall Width: 640mm * Upholstery Nylon (Flame Retardant) * Colour: Silver Warranty: 12 months on frame & crossbars



Yes

Self-Propelled 20" LINK
SPECIAL ORDER

Days Premium Link Wheelchair is a lightweight adaptable wheelchair suitable for full-time users and is designed for maximum practicality and adjustability.

Light weight aluminium - 17.2kg

- Quick release wheels
- Solid tyres
- 16", 18", 20" and 22" seat widths available
- Swing in/out legrests
- Swing back removable armrests
- Height adjustable back upholstery
- Adjustable castors
- Multi position rear wheels
- Drum brakes on the 20"-22"
- Anti-tipper comes as standard inclusion





Yes

Self-Propelled Bariatric 20"
MERITS

Reinforced double cross-brace for extra strength with weight capacity of 182kg
Detachable full length armrests
Detachable swing away leg rests with oversize footplates, with optional elevating leg rests
24" x 1.3/8" puncture free foam filled tyres

Overall Width 73.5cm(20")
Overall Length 123cm
Overall Height 91 cm
Seat Width 48cm(20")
Seat Depth 46cm
Seat height from ground 49cm
Seat Back Height 40cm
Footrest Length
Adjustment 37cm – 50cm
Armrest Height 23cm
Weight Capacity 182kg
Total weight 15.5kg(20")



Yes	Self-Propelled Bariatric 22" MERITS	<p>Reinforced double cross-brace for extra strength with weight capacity of 182kg</p> <p>Detachable full length armrests</p> <p>Detachable swing away leg rests with oversize footplates, with optional elevating leg rests</p> <p>24" x 1.3/8" puncture free foam filled tyres</p>	<p>Overall Width 78.5cm(22")</p> <p>Overall Length 123cm</p> <p>Overall Height 91 cm</p> <p>Seat Width 56cm(22")</p> <p>Seat Depth 46cm</p> <p>Seat height from ground 49cm</p> <p>Seat Back Height 40cm</p> <p>Footrest Length Adjustment 37cm – 50cm</p> <p>Armrest Height 23cm</p> <p>Weight Capacity 182kg</p> <p>Total weight 16.5kg(22")</p>	
Yes	Self-Propelled Breezy 16" with solid Tyres and Fixed armrests.	<p>The Breezy BasiX² is a lightweight folding wheelchair that features amazing adjustability and rigidity, combined with a modern and clean look. Adjustments and stump supports available for amputee set up.</p> <p>Recent Improvements include greater seat depth adjustment (41cm to 46cm), a tension adjustable backrest which is extendable from 41cm to 46cm and a more robust, rigid castor attachment (0, 3 and 6 degrees).</p>	<p>Seat Width: 38 - 52cm</p> <p>Seat Depth: 41 - 46cm</p> <p>Seat Height: 42.5 - 51cm</p> <p>Seat Angle: 0 - 6 degrees</p> <p>Backrest Height: 41 - 46cm</p> <p>SWL: 125kg</p> <p>Total Weight: 14.9Kg</p>	

Yes

Self-Propelled Breezy 16" with pneumatic Tyres and Fixed armrests.

The Breezy BasiX² is a lightweight folding wheelchair that features amazing adjustability and rigidity, combined with a modern and clean look. Adjustments and stump supports available for amputee set up.

Recent Improvements include greater seat depth adjustment (41cm to 46cm), a tension adjustable backrest which is extendable from 41cm to 46cm and a more robust, rigid castor attachment (0, 3 and 6 degrees).

Seat Width: 38 - 52cm
Seat Depth: 41 - 46cm
Seat Height: 42.5 - 51cm
Seat Angle: 0 - 6 degrees
Backrest Height: 41 - 46cm
SWL: 125kg
Total Weight: 14.9Kg



Yes

Self-Propelled Breezy 16" with solid Tyres and Adjustable armrests.

The Breezy BasiX² is a lightweight folding wheelchair that features amazing adjustability and rigidity, combined with a modern and clean look. Adjustments and stump supports available for amputee set up.

Recent Improvements include greater seat depth adjustment (41cm to 46cm), a tension adjustable backrest which is extendable from 41cm to 46cm and a more robust, rigid castor attachment (0, 3 and 6 degrees).

Seat Width: 38 - 52cm
Seat Depth: 41 - 46cm
Seat Height: 42.5 - 51cm
Seat Angle: 0 - 6 degrees
Backrest Height: 41 - 46cm
SWL: 125kg
Total Weight: 14.9Kg



Yes

Self-Propelled Breezy 18" with solid Tyres and Fixed armrests.

The Breezy BasiX² is a lightweight folding wheelchair that features amazing adjustability and rigidity, combined with a modern and clean look. Adjustments and stump supports available for amputee set up.

Recent Improvements include greater seat depth adjustment (41cm to 46cm), a tension adjustable backrest which is extendable from 41cm to 46cm and a more robust, rigid castor attachment (0, 3 and 6 degrees).

Seat Width: 38 - 52cm
Seat Depth: 41 - 46cm
Seat Height: 42.5 - 51cm
Seat Angle: 0 - 6 degrees
Backrest Height: 41 - 46cm
SWL: 125kg
Total Weight: 14.9Kg



Yes

Self-Propelled Breezy 18" with pneumatic Tyres and Fixed armrests.

The Breezy BasiX² is a lightweight folding wheelchair that features amazing adjustability and rigidity, combined with a modern and clean look. Adjustments and stump supports available for amputee set up.

Recent Improvements include greater seat depth adjustment (41cm to 46cm), a tension adjustable backrest which is extendable from 41cm to 46cm and a more robust, rigid castor attachment (0, 3 and 6 degrees).

Seat Width: 38 - 52cm
Seat Depth: 41 - 46cm
Seat Height: 42.5 - 51cm
Seat Angle: 0 - 6 degrees
Backrest Height: 41 - 46cm
SWL: 125kg
Total Weight: 14.9Kg



Yes

Self-Propelled Breezy 18" with solid Tyres and Adjustable armrests.

The Breezy BasiX² is a lightweight folding wheelchair that features amazing adjustability and rigidity, combined with a modern and clean look. Adjustments and stump supports available for amputee set up.

Recent Improvements include greater seat depth adjustment (41cm to 46cm), a tension adjustable backrest which is extendable from 41cm to 46cm and a more robust, rigid castor attachment (0, 3 and 6 degrees).

SWL: 125kg
Weight: 14.9kg
Seat Width: 38 - 52cm
Seat Height: 42.5 - 51cm
Seat Depth: 41 - 46cm



Transit Wheelchairs

Yes

Transit Wheelchair 14"
ASPIRE

The Transit 2 – is designed for an individual or facility who requires a robust chair that will meet the rigorous demands of the community and healthcare environments.
Swing away and removable arm rests – this facilitates user transfers and also makes for easier transport and storage.

SWL: 140kg
Weight: 16.4kg



Yes

Transit Wheelchair 16"
ASPIRE

The Transit 2 – is designed for an individual or facility who requires a robust chair that will meet the rigorous demands of the community and healthcare environments.
Swing away and removable arm rests – this facilitates user transfers and also makes for easier transport and storage.

SWL: 140kg
Weight: 16.6kg



Yes

Transit Wheelchair 16" ESCAPE

With a foldable aluminium frame, the Escape has been specially designed for ease of transportation and storage, making it an ideal chair for occasional use. Suitable for both indoors and out, it gives you flexibility wherever you go.

In addition to being portable, these attendant-propelled chairs offer both comfort and practicality. With a full-length armrest, padded upholstery and detachable swing-away footrests, they provide all the basic comfort necessities. Peace of mind and low maintenance are provided by the puncture-proof tyres.

Seat Height - Front (from the ground): 510mm
Seat Height - Rear (from the ground): 490
Seat Width (on the outside of the tubes): 450mm
Seat Width (on the inside of the tubes): 400mm
Seat Depth: 400mm
Width Between Arms: 400mm
Armrest Height Above Seat: 230mm
Backrest Height (from the seat): 400mm
Footrest Length (from seat to surface of footrest): 350-450mm





Yes

Transit Wheelchair 18inch
SWIFT

The Days Swift Transit Attendant Propelled Wheelchair has been designed to be a light, practical and easy to transport wheelchair, making it an ideal occasional use chair for both indoor and outdoor use. Its compact dimensions allow good manoeuvrability indoors & the light construction and ease to fold for storage or transport make it a very versatile chair. Attendant propelled in one size.

Light frame 7.8kg
Folding backrest (for storage or transport)
Fixed arms
Swing away and removable leg rests
Solid tyres (no maintenance required)
Extended brake handle (ease of brake application)



Yes	Transit Wheelchair 18" WHIRL	<p>The Healthcare Whirl is a versatile wheelchair with a range of features to accommodate different user needs. * This sleek and stylish wheelchair is practical and will look good when you are out and about. Features: * Quick release wheels (for transport) * Removable armrests * Swing away and removable leg rests * Solid tyres (no maintenance required)</p>	<p>Specifications: * Weight: 16.8kg (complete) * User Weight: 120kg * Seat Width: 450mm * Seat Depth: 430mm * Backrest Height: 440mm * Seat Height - Front: 500mm * Seat Height - Rear: 475mm * Front Castor Wheel: 8" x 1" PU * Rear Wheel: 24" x 1 3/8" PU * Overall Width: 640mm * Upholstery Nylon (Flame Retardant) * Colour: Silver Warranty: 12 months on frame & crossbars</p>	
Yes	Transit Wheelchair 20" ASPIRE	<p>The Transit 2 – is designed for an individual or facility who requires a robust chair that will meet the rigorous demands of the community and healthcare environments. Swing away and removable arm rests – this facilitates user transfers and also makes for easier transport and storage.</p>	<p>SWL: 140kg Weight: 17.9kg</p>	

Stump Supports

Yes

Yes

BREEZY Stump Support -
RIGHT




This wheelchair Accessory
stump support replaces the
footrest in the Breezy Basix
wheelchair.





BREEZY Stump Support -
LEFT

This wheelchair Accessory
stump support replaces the
footrest in the Breezy Basix
wheelchair.



CHSP EQUIP IN SCOPE LIST	PRODUCT	PRODUCT DESCRIPTION	SPECIFICATIONS	IMAGE
Yes	CHAIRS AND SEATING			
	Foot Stools			
	Foot Stool / Leg Rest	The Legrest Stool is height adjustable, designed with a 10 degree angle for comfort and support and should be placed with the higher end against the chair.	Seat Pad 400mm x 500mm x 40mm Height Adjustable 420mm - 570mm Maximum User Weight 160kg	
	Kitchen Stools			
	Kitchen Perching Stool - With Arms	Ideal for use in the kitchen, or any activity where a perched seating position is of assistance Height adjustable Padded seat and backrest Integrated arm supports assist transfers Vinyl upholstery	Width Between Arms: 455mm Seat Depth: 300mm Seat to Top of Backrest: 300mm Leg Height Adjustment: 500mm Floor to Seat: 635mm Maximum User Weight: Up to 100kg	
	Utility Chairs			
	Low Back Utility Chair	The Standard Utility Chair features a contoured backrest, padded armrests, vinyl upholstery, adjustable height, and a powder coated steel frame. The adjustable height assists users who struggle to get out of standard chairs.	Overall Width 580mm Seat Width 455mm Seat Depth 460mm Seat Height 450 - 600mm Back Height 465mm Maximum User Weight 160kg	

Yes	High Back Utility Chair	This chair can be used in a lounge room instead of a couch or regular armchair that is too low to get in and out of with ease. Ffeatures a contoured backrest, padded armrests, vinyl upholstery, adjustable height, and a powder coated steel frame. Available in Teal/Blue and Mocha.	Overall Width 640mm Seat Width 490mm Seat Depth 480mm Seat Height 455 - 605mm Back Height 730mm Maximum User Weight 160kg	
	High Back Utility Chair - Bariatric	This chair can be used in a lounge room instead of a couch or regular armchair that is too low to get in and out of with ease. Wider seat widths and reinforced frames to suit the bariatric client Adjustable height and seat depth Contoured backrest Powder coated steel frame	Upholstery Champagne Vinyl Seat Width 600mm Seat Depth 490-550mm Seat Height 430-580mm Arm Height 250mm Back Height 740mm Overall Width 800mm Max User Weight 300kg	
	Adult Postural Chairs			

No

Chair Air Deluxe Bed Adult with
Footrest - Large

Cushioned adjustable air cell support has been incorporated into the design. The four lateral support wings fold away for easy transferring of the patient. This chair has a reclinable seat along with an adjustable backrest and calf rest, tray, neck support and footplate. The frame is completely constructed from sturdy stainless steel.



Features:

- Adjustable calf rest
- Reclinable seat
- Air system for pressure relief
- Full stainless steel frame
- Three lockable castors and one steering castor (non-marking)
- Breathable polyurethane cover (fire retardant)
- Replaceable and easy to clean cover
- Four folding wings
- Standard dining tray
- Magazine pocket

Seat Height: 63cm
Seat Width: 51cm
Seat Depth: 50cm
Overall Width: 73cm
Overall Length: 135 - 188cm
Table to Floor: 88cm
Table Slide: 23cm
Table Size: 60 x 40.5 x 2cm
Back Reclining: 102° - 136°
Seat Tilting: 0 - 18°
Leg Elevating: 70° down to 5° up
Weight of Chair: 43.5kg
Maximum User Weight: 120kg
Warranty: 1 year



Electric Recliner Chairs

No	Monarch Chair (PETITE SINGLE MOTOR)	<p>Seamed backrest provides additional lumbar support</p> <p>Padded legrest is synchronised with the backrest</p> <p>Full chaise pad provides additional legrest support and comfort</p> <p>Supplied with magazine pocket</p> <p>Hand held controller is simple to use</p> <p>Rise or recline motion stops at any point when button is released</p>	<p>Petite</p> <p>Seat Width 480mm</p> <p>Seat Depth 400mm</p> <p>Seat Height 460mm</p> <p>Arm Height 140mm</p> <p>Back Height 620mm</p> <p>Overall Width 780mm</p> <p>Max User Weight 130Kg</p>	
No	Monarch Chair (SINGLE MOTOR)	<p>Seamed backrest provides additional lumbar support</p> <p>Padded legrest is synchronised with the backrest</p> <p>Full chaise pad provides additional legrest support and comfort</p> <p>Supplied with magazine pocket</p> <p>Hand held controller is simple to use</p> <p>Rise or recline motion stops at any point when button is released</p>	<p>Standard</p> <p>Seat Width 550mm</p> <p>Seat Depth 460mm</p> <p>Seat Height 500mm</p> <p>Arm Height 190mm</p> <p>Back Height 680mm</p> <p>Overall Width 820mm</p> <p>Max User Weight 130Kg</p>	

No

Monarch Chair (DUAL MOTOR)

Seamed backrest provides additional lumbar support
Padded legrest is synchronised with the backrest
Full chaise pad provides additional legrest support and comfort
Supplied with magazine pocket
Hand held controller is simple to use
The dual motor provides independent control of backrest and legrest.

Standard
Seat Width 550mm
Seat Depth 460mm
Seat Height 500mm
Arm Height 190mm
Back Height 680mm
Overall Width 820mm
Max User Weight 130Kg



No

Mayfair Chair - Black (SINGLE
MOTOR - SPACE SAVER)

The Select Mayfair Chair features a very comfortable 3 Section Waterfall Backrest Design for Optimum Comfort and Upper, Middle and Lower back support.

Quality Foams and no Sag springs provide supreme Postural Support.

Backup Power Supply.

Featuring quality Actuators and American Branded Mechanisms means quiet motors for a Smooth & Discreet Operation.

Standard
Seat Width 520mm
Seat Depth 480mm
Seat Height 480mm
Overall Width 850mm
Max User Weight 150Kg



No

Mayfair Chair - Vanilla (SINGLE
MOTOR - SPACE SAVER)

The Select Mayfair Chair features a very comfortable 3 Section Waterfall Backrest Design for Optimum Comfort and Upper, Middle and Lower back support.

Quality Foams and no Sag springs provide supreme Postural Support.

Backup Power Supply.

Featuring quality Actuators and American Branded Mechanisms means quiet motors for a Smooth & Discreet Operation.

Standard
Seat Width 520mm
Seat Depth 480mm
Seat Height 480mm
Overall Width 850mm
Max User Weight 150Kg



No

Drive Portland Recliner PETITE
(DUAL MOTOR)

The Portland Riser Recliner Lift Chair is designed with comfort in mind. It features an exclusive memory foam seat and its dual motor mechanism allows the user to find their perfect comfort position in minutes.

Height: 101cm
Width: 83cm
Depth: 90cm
Seat Height: 43cm
Seat Width: 47cm
Seat Depth: 48cm
SWL: 150kg

The Portland Riser Lift Chair Armchair features luxurious Sofia fabric and is available in two colours- Mushroom and Graphite.

Features of the Portland Riser Lift Chair include:

Reclines to an excellent snooze position.
High quality upholstery.
Excellent quiet and smooth actions.
Built in convenient storage pockets.
Battery back up facility
New 2 year warranty
Weight capacity 150kg



No

Drive Portland Recliner (DUAL MOTOR)

The Portland Riser Recliner Lift Chair is designed with comfort in mind. It features an exclusive memory foam seat and its dual motor mechanism allows the user to find their perfect comfort position in minutes.

Height: 110cm
Width: 86cm
Depth: 94cm
Seat Height: 438cm
Seat Width: 51cm
Seat Depth: 54cm
SWL: 150kg

The Portland Riser Lift Chair Armchair features luxurious Sofia fabric and is available in two colours- Mushroom and Graphite.

Features of the Portland Riser Lift Chair include:

Reclines to an excellent snooze position.
High quality upholstery.
Excellent quiet and smooth actions.
Built in convenient storage pockets.
Battery back up facility
New 2 year warranty
Weight capacity 150kg



Step Stools



No





Step Stool with optional hand rail.

Chrome plated steel frame step stool with hand rail. It has a slip-resistant fluted rubber covered step and a tall hand rail can be fitted to one side. It also has non-marking, slip resistant rubber tips.

Height of Step: 23cm
• Height of Handrail (from floor): 104cm
• Overall Width: 44cm
• Overall Depth: 35.5cm
• Step Width: 36cm
• Step Depth: 29cm
• Max User Weight: 160kg



CHSP EQUIP IN SCOPE LIST	PRODUCT	PRODUCT DESCRIPTION	SPECIFICATIONS	IMAGE
Yes	CONTINENCE			
	Bed Pans			
	Bed Pan - Slipper	Assists toileting of user whilst confined to bed Wedge shape assists in correct positioning Poly Propylene construction		
Yes	Bed Pan Plastic	Light weight polypropylene bed pan for toileting while remaining in bed Wide opening cut-out area to allow for correct positioning Durable moulded plastic design with autoclavable construction Compatible with most Mobile Shower Commodes		
Urinals				

Yes	Urinal Male - Non spill	<p>Autoclavable for sanitising</p> <p>Integrated handle</p> <p>Patented non-spill spout prevents spillage in any position</p> <p>Spout unscrews easily for disposal and cleaning</p> <p>100ml graduation markings</p> <p>Moulded plastic</p>	<p>Weight 240g</p> <p>1 litre capacity</p> <p>500ml non-spill capacity</p>	
Yes	Urinal Male - Standard	<p>Standard urinal bottle for males. Ideal for bed confinement or travelling</p>	<p>Graduated up to 1000 ml</p> <p>Secured lid to reduce spills and odours</p>	
Yes	Urinal Female	<p>Assists the user or carer in personal toileting</p> <p>Ideal for bed confinement or travelling</p> <p>Anatomically shaped for females</p> <p>Secured lid to reduce spills and odours</p>	<p>Graduated up to 1000 ml</p>	
Yes	<p>Absorbent Sheets</p> <p>Absorbent Sheet - With Bed Tuck Ins</p>	<p>Highly absorbent draw sheet with side tuck-ins. This item comes with built in waterproofing.</p>	<p>Width (not incl. tuck-ins): 1000mm</p> <p>Length: 1000mm</p> <p>Absorption Capacity: 2.5L</p>	

Yes

McIntoshes

Waterproof Sheeting / McIntosh Mack-E Mattress Protector with wings
PVC



CHSP EQUIP IN
SCOPE LIST

PRODUCT

PRODUCT DESCRIPTION

SPECIFICATIONS

IMAGE

ELECTRIC BEDS**Electric Adjustable Beds**

No

Community Electric Bed

A versatile and adaptable nursing bed which is easy to use and built for strength and durability for the home care market. Electric bed, with electric backrest/headraise, knee break and manual leg lift. Trendelenburg and Reverse Trendelenburg. 4 locking castors. Bed folds for transportation.

- Easy to fold and superior clearance when folded for transport or storage
- Multiple positioning to suit a range of care needs
- Easy to use hand controller with the ability to lockout individual functions
- Linak actuators providing superior reliability
- Strong and durable frame with a 5-year warranty

- Safe Working Load - 220kg
- Dimensions - 2324x1050mm
- Height Range - 330 to 730mm
- Mattress size - 1960
- Bed Frame Weight - 74 kg
- Warranty - 5 year bed frame / 2 year electrical / 1 year castor



No

Electric Bed, Aged Care

A fully-featured, easy to operate, electric bed.
Four-part, high low profiling bed
Adjustable bed height, backrest, knee break and heel lift
Quick-release head and foot panels
Choice of panel designs
Central locking castors
In-built 200mm base extension
Electric hand controller with auto contour and a lockout function
4 Detachable mattress retainers
Drainage bag holder

- Safe working load - 220 kg
- Maximum patient weight - 185 kg
- Warranty - 3 year Electrical / 10 year Frame / 1 year Back-Up Battery
- Frame material - Powder Coated Steel
- Model - Standard, King
- Overall dimensions
- W 920 mm† L 2036 mm‡
- W 1090 mm† L 2036 mm‡
- Height base - 330 - 780 mm 330 - 780 mm
- Mattress dimensions
recommended Single - W 860 - 900 mm
- Mattress dimensions
recommended Single - L 1960 - 2200* mm
- Mattress dimensions
recommended King - W 1050 - 1070 mm



No

iCare King Single Complete Bed
Package ONYX

Safe working load 200kg on king single.
Hand Controller Backlit for night time use. Easy press buttons. Wired stretch cord. Hook attachment for hanging to side of bed. Wipeable.
Hard wearing commercial grade fabric. Onyx range (BLACK) 82% Polyester, 18% Olefin. Stone range (GREY) 70% Polyester, 30% Acrylic.
Strong metal mesh under mattress allowing increased airflow
8x heavy duty swivel castors with one lockable castor on each corner.
Head lift angle 62° and knee break angle 42°
Max height 66cm from floor to top of base
Min height 22cm from floor to top of base.
Includes the IC20 Mattress a mid range mattress with a soft feel but excellent support and a matching headboard in ONYX fabric.

Safe working load -
200kg.
Minimum height from
floor to base is 220mm.



No

iCare Long Single Complete
Bed Package STONE

Safe working load 200kg on long single.
Hand Controller Backlit for night time use. Easy press buttons. Wired stretch cord. Hook attachment for hanging to side of bed. Wipeable.
Hard wearing commercial grade fabric. Onyx range (BLACK) 82% Polyester, 18% Olefin. Stone range (GREY) 70% Polyester, 30% Acrylic.
Strong metal mesh under mattress allowing increased airflow
8x heavy duty swivel castors with one lockable castor on each corner.
Head lift angle 62° and knee break angle 42°
Max height 66cm from floor to top of base
Min height 22cm from floor to top of base.
Includes the IC20 Mattress a mid range mattress with a soft feel but excellent support and a matching headboard in STONE fabric.

Safe working load -
200kg.
Minimum height from
floor to base is 220mm.



No

iCare King Single Complete Bed
Package STONE

Safe working load 200kg on king single.
Hand Controller Backlit for night time use. Easy press buttons. Wired stretch cord. Hook attachment for hanging to side of bed. Wipeable.
Hard wearing commercial grade fabric. Onyx range (BLACK) 82% Polyester, 18% Olefin. Stone range (GREY) 70% Polyester, 30% Acrylic.
Strong metal mesh under mattress allowing increased airflow
8x heavy duty swivel castors with one lockable castor on each corner.
Head lift angle 62° and knee break angle 42°
Max height 66cm from floor to top of base
Min height 22cm from floor to top of base.
Includes the IC20 Mattress a mid range mattress with a soft feel but excellent support and a matching headboard in STONE fabric.

Safe working load -
200kg.
Minimum height from
floor to base is 220mm.



Accessories

No

Bedstick Low - Townsend for Electric Beds

Clamped firmly to bed base railing
Acts as a support to roll over, to sit from lying position and to assist with getting in and out of bed.

For use for Hospital Style Beds

Manufactured in Adelaide,
Powdercoated mild steel. Clean using mild detergent and water and dry with soft cloth.

Low profile 525 x 175
2.1kg
Max User Weight 136kg



No


BedstickStandard - Townsend for Electric Beds

Clamped firmly to bed base railing
Acts as a support to roll over, to sit from lying position and to assist with getting in and out of bed.

For use for Hospital Style Beds

Manufactured in Adelaide,
Powdercoated mild steel. Clean using mild detergent and water and dry with soft cloth.

Standard Height 715x175
2.4kg
Max User Weight 136kg

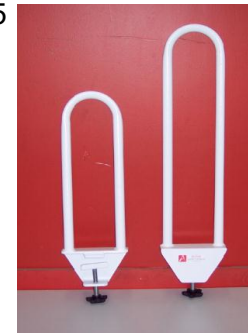

Discuss with CL

NEW Lawson Bedstick with
11cm gap between rails

Installed between mattress and base
Acts as a support roll over, to sit from a
lying position and to assist with getting
in and out of bed
For use to suit domestic ensemble or
standard mattress beds

NEW 11cm ap between



No

Bedrail for Elec Bed - Fold
Down - 3/4 length

Clamp on, fold down metal side rail
system
attaches firmly to bed frame. Three
quarter
length rails for increased access.



No

Bed Rail - Full Length for
Electric Bed

Drop down metal side rail system
attaches firmly
to bed frame. Full length rails for added
safety
and security.



No

Bed Rail Cover - Padded for Electric Bed

Soft cushions to encase metal bed rails Standard 670 x 1260mm protects limbs from hitting rails or getting caught.
Soft, medical grade foam cushions, stretch cover, easy secure Velcro closure, easy wipe clean.
Two sizes available, fits majority of bed rails.



No

Pole Self Help For Community Bed

Support pole to aid patient repositioning while lying in bed.
Mounts into inbuilt brackets at either side of bed head.



CHSP EQUIP
IN SCOPE LIST

PRODUCT

PRODUCT DESCRIPTION

SPECIFICATIONS

IMAGE

No

PRESSURE CARE

Pressure Cushions - Air

Cushion ROHO Low Profile 8x9
(fits 14x16 W/C)

The Low-Profile Cushions effectively redistribute pressure and conform to the user's body in order to prevent against ulcers and tissue breakdown. Its soft, air-filled cells permit blood flow while simultaneously providing users with the additional support they require while seated. The ROHO Low-Profile Cushion maximizes stability while reducing the risk of potential tissue breakdown.

Available with either a nylon lycra or incontinent resistant cover. Both are machine washable, and they protect the cushion's surface from wear and tear as well as provide aid during slide transfers. The bottom of the cover has a non-skid coating to prevent sliding while in use. The sturdy neoprene rubber cells resist liquid/moisture and clean easily with soap and water.

Cell size 8 x 9
Dimensions 38 x 43 cm
Profile Low
Height 60mm
Max User Weight
Unlimited



No

Cushion ROHO Low Profile 9x9

The Low-Profile Cushions effectively redistribute pressure and conform to the user's body in order to prevent against ulcers and tissue breakdown. Its soft, air-filled cells permit blood flow while simultaneously providing users with the additional support they require while seated. The ROHO Low-Profile Cushion maximizes stability while reducing the risk of potential tissue breakdown.

Available with either a nylon lycra or incontinent resistant cover. Both are machine washable, and they protect the cushion's surface from wear and tear as well as provide aid during slide transfers. The bottom of the cover has a non-skid coating to prevent sliding while in use. The sturdy neoprene rubber cells resist liquid/moisture and clean easily with soap and water.

Cell size 9 x 9
Dimensions 43 x 43 cm
Profile Low
Height 60mm
Max User Weight
Unlimited



No

Cushion ROHO Low Profile 10 cells x 9 cells (fits 18x16 W/C)

The Low-Profile Cushions effectively redistribute pressure and conform to the user's body in order to prevent against ulcers and tissue breakdown. Its soft, air-filled cells permit blood flow while simultaneously providing users with the additional support they require while seated. The ROHO Low-Profile Cushion maximizes stability while reducing the risk of potential tissue breakdown.

Available with either a nylon lycra or incontinent resistant cover. Both are machine washable, and they protect the cushion's surface from wear and tear as well as provide aid during slide transfers. The bottom of the cover has a non-skid coating to prevent sliding while in use. The sturdy neoprene rubber cells resist liquid/moisture and clean easily with soap and water.

Cell size 10 x 9
Dimensions 47 x 43cm
Profile Low
Height 60mm
Weight Unlimited



No

Cushion ROHO Low Profile 10
cells x 10 cells

The Low-Profile Cushions effectively redistribute pressure and conform to the user's body in order to prevent against ulcers and tissue breakdown. Its soft, air-filled cells permit blood flow while simultaneously providing users with the additional support they require while seated. The ROHO Low-Profile Cushion maximizes stability while reducing the risk of potential tissue breakdown.

Available with either a nylon lycra or incontinent resistant cover. Both are machine washable, and they protect the cushion's surface from wear and tear as well as provide aid during slide transfers. The bottom of the cover has a non-skid coating to prevent sliding while in use. The sturdy neoprene rubber cells resist liquid/moisture and clean easily with soap and water.

Cell size 10 x 10
Dimensions 47 x 47 cm
Profile Low
Height 60mm
Max User Weight
Unlimited



No

Cushion ROHO Low Profile 11 cells x 10 cells (fits 20x18 W/C)

The Low-Profile Cushions effectively redistribute pressure and conform to the user's body in order to prevent against ulcers and tissue breakdown. Its soft, air-filled cells permit blood flow while simultaneously providing users with the additional support they require while seated. The ROHO Low-Profile Cushion maximizes stability while reducing the risk of potential tissue breakdown.

Available with either a nylon lycra or incontinent resistant cover. Both are machine washable, and they protect the cushion's surface from wear and tear as well as provide aid during slide transfers. The bottom of the cover has a non-skid coating to prevent sliding while in use. The sturdy neoprene rubber cells resist liquid/moisture and clean easily with soap and water.

Cell size 11 x 10
Dimensions 52 x 47 cm
Profile Low
Height 60mm
Max User Weight Unlimited



No

Cushion ROHO High Profile 10
cells x 10 cells

The High Profile Cushions redistribute pressure and conform to the user's body in order to prevent against ulcers and tissue breakdown. Its soft, air-filled cells permit blood flow while simultaneously providing users with the additional support they require while seated. The ROHO High-Profile Cushion assists those who are unable to independently shift their own weight.

Available with your choice of either a nylon lycra or incontinent resistant cover. Both are machine washable, and they protect the cushion's surface from wear and tear as well as provide aid during slide transfers. The bottom of the cover has a non-skid coating to prevent sliding while in use.

Easy Clean Wheelchair Cushion

The sturdy neoprene rubber cells resist liquid/moisture and clean easily with soap and water.

Cell size 10 x 10

Dimensions 47 x 47 cm

Profile High

Height 100mm

Max User Weight

Unlimited



No

Cushion ROHO High Profile 11 cells x 10 cells (fits 20x18 W/C)

The High Profile Cushions redistribute pressure and conform to the user's body in order to prevent against ulcers and tissue breakdown. Its soft, air-filled cells permit blood flow while simultaneously providing users with the additional support they require while seated. The ROHO High-Profile Cushion assists those who are unable to independently shift their own weight.

Available with your choice of either a nylon lycra or incontinent resistant cover. Both are machine washable, and they protect the cushion's surface from wear and tear as well as provide aid during slide transfers. The bottom of the cover has a non-skid coating to prevent sliding while in use.

Easy Clean Wheelchair Cushion

The sturdy neoprene rubber cells resist liquid/moisture and clean easily with soap and water.

Cell size 11 x 10
Dimensions 52 x 47 cm
Profile High
Height 100mm
Max User Weight Unlimited
Screen reader support enabled.



No

Cushion ROHO High Profile 8x9

The High Profile Cushions redistribute pressure and conform to the user's body in order to prevent against ulcers and tissue breakdown. Its soft, air-filled cells permit blood flow while simultaneously providing users with the additional support they require while seated. The ROHO High-Profile Cushion assists those who are unable to independently shift their own weight.

Available with your choice of either a nylon lycra or incontinent resistant cover. Both are machine washable, and they protect the cushion's surface from wear and tear as well as provide aid during slide transfers. The bottom of the cover has a non-skid coating to prevent sliding while in use.

Easy Clean Wheelchair Cushion

The sturdy neoprene rubber cells resist liquid/moisture and clean easily with soap and water.

Cell size 8 x 9

Dimensions 38 x 43 cm

Profile High

Height 100mm

Max User Weight

Unlimited



No

Cushion ROHO High Profile 9x9

The High Profile Cushions redistribute pressure and conform to the user's body in order to prevent against ulcers and tissue breakdown. Its soft, air-filled cells permit blood flow while simultaneously providing users with the additional support they require while seated. The ROHO High-Profile Cushion assists those who are unable to independently shift their own weight.

Available with your choice of either a nylon lycra or incontinent resistant cover. Both are machine washable, and they protect the cushion's surface from wear and tear as well as provide aid during slide transfers. The bottom of the cover has a non-skid coating to prevent sliding while in use.

Easy Clean Wheelchair Cushion

The sturdy neoprene rubber cells resist liquid/moisture and clean easily with soap and water.

Cell size 9 x 9

Dimensions 43 x 43 cm

Profile High Height

100mm Max User

Weight Unlimited



No

Cushion ROHO Quadtro 8 cells
x 9 cells HP (fits 14x16 W/C)

The Quadtro Select cushion is made for individuals who require superior pressure relief and protection from soft tissue breakdown. Its high-profile nature facilitates blood flow via advanced shape-conforming technology.

The ROHO Quadtro Select is lightweight, easy to transport, and easy to clean (just use soap and water). It features four individual sections that are uniformly inflated through a single valve. A simple push of the ISOFLO Memory Control knob locks air into each section, therefore filling them to maximum capacity. The ROHO cushion is also available in varying dimensions to accommodate for users' individual seating needs.

Cushion Cover Included

This seat cushion comes with a durable cover that reduces unnecessary friction and shear. Spill something on your seat? Machine wash the cover in cold water for an easy clean up.

Cell size 8 x 9
Dimensions 38 x 43 cm
Profile High
Height 100mm
Max User Weight
Unlimited



No

Cushion ROHO HP Quadtro 9 cells x 9 cells (fits 16x16 W/C)

The Quadtro Select cushion is made for individuals who require superior pressure relief and protection from soft tissue breakdown. Its high-profile nature facilitates blood flow via advanced shape-conforming technology.

The ROHO Quadtro Select is lightweight, easy to transport, and easy to clean (just use soap and water). It features four individual sections that are uniformly inflated through a single valve. A simple push of the ISOFLO Memory Control knob locks air into each section, therefore filling them to maximum capacity. The ROHO cushion is also available in varying dimensions to accommodate for users' individual seating needs.

Cushion Cover Included

This seat cushion comes with a durable cover that reduces unnecessary friction and shear. Spill something on your seat? Machine wash the cover in cold water for an easy clean up.

Cell size 9 x 9
Dimensions 43 x 43 cm
Profile High
Height 100mm
Max User Weight
Unlimited



No

Cushion ROHO Quadtro 10 cells x 9 cells High Profile (fits 18x16 W/C)

The Quadtro Select cushion is made for individuals who require superior pressure relief and protection from soft tissue breakdown. Its high-profile nature facilitates blood flow via advanced shape-conforming technology.

The ROHO Quadtro Select is lightweight, easy to transport, and easy to clean (just use soap and water). It features four individual sections that are uniformly inflated through a single valve. A simple push of the ISOFLO Memory Control knob locks air into each section, therefore filling them to maximum capacity. The ROHO cushion is also available in varying dimensions to accommodate for users' individual seating needs.

Cushion Cover Included

This seat cushion comes with a durable cover that reduces unnecessary friction and shear. Spill something on your seat? Machine wash the cover in cold water for an easy clean up.

Cell size 10 x 9
Dimensions 47 x 43 cm
Profile High
Height 100mm
Max User Weight
Unlimited



No

Cushion ROHO HP Quadtro
10x10

The Quadtro Select cushion is made for individuals who require superior pressure relief and protection from soft tissue breakdown. Its high-profile nature facilitates blood flow via advanced shape-conforming technology.

The ROHO Quadtro Select is lightweight, easy to transport, and easy to clean (just use soap and water). It features four individual sections that are uniformly inflated through a single valve. A simple push of the ISOFLO Memory Control knob locks air into each section, therefore filling them to maximum capacity. The ROHO cushion is also available in varying dimensions to accommodate for users' individual seating needs.

Cushion Cover Included

This seat cushion comes with a durable cover that reduces unnecessary friction and shear. Spill something on your seat? Machine wash the cover in cold water for an easy clean up.

Cell size 10 x 10
Dimensions 47 x 47 cm
Profile High
Height 100mm
Max User Weight
Unlimited



No

Cushion HP ROHO Quadtro
11x10 (fits 20x18 W/C)

The Quadtro Select cushion is made for individuals who require superior pressure relief and protection from soft tissue breakdown. Its high-profile nature facilitates blood flow via advanced shape-conforming technology.

The ROHO Quadtro Select is lightweight, easy to transport, and easy to clean (just use soap and water). It features four individual sections that are uniformly inflated through a single valve. A simple push of the ISOFLO Memory Control knob locks air into each section, therefore filling them to maximum capacity. The ROHO cushion is also available in varying dimensions to accommodate for users' individual seating needs.

Cushion Cover Included

This seat cushion comes with a durable cover that reduces unnecessary friction and shear. Spill something on your seat? Machine wash the cover in cold water for an easy clean up.

Cell size 11 x 10
Dimensions 52 x 47 cm
Profile High
Height 100mm
Max User Weight
Unlimited



Combination Cushions

No

Cushion Jay 2 Standard
460mmx410mm (18"x16")

Contoured foam cushion with Jay Flow 46 x 41 cm
fluid pad 18"x16"
Designed for the client with moderate
to high postural and pressure care
needs
Firm, anatomically correct contoured
foam base creates a foundation for
postural stability
Provides excellent lateral and
forward/rear stability
Fluid Tripad in rear seat well maintains
even pressure distribution and keeps
fluid under areas at high risk of skin
breakdown
Supplied with breathable ballistic
stretch cover



No

Cushion ROHO Hybrid Elite
460mmx460mm (fits 18x16
W/C)

Combines ROHOs Air Flotation
qualities and Jays base technology
Ideal for clients who require moderate
to high pressure care and positional
support needs
Contoured foam base improves
alignment of the pelvis and lower
extremities, enhancing sitting posture
and function
Air cells are interconnected and their
unique shape promote blood circulation
to facilitate pressure prevention,
healing and relief
Flexible and adjustable air cells adapt
to patient contours and reduce shear
forces while transferring or
repositioning
Durable neoprene construction is soft
to suit skin sensitivity, and is easily
cleaned
Jay foam base is made from fluid-
resistant closed cell foam for additional
stability, assisting in safer wheelchair
transfers
Two-way stretch cover with encased
comfort foam is machine washable
Supplied with pump
Cushions is available with fluid resistant
cover

Dimensions 46 x 46cm
Height 95mm
Max User Weight 225kg



Foam Cushions

Yes

Cushion Basic Foam 405mm x
405mm Black Cover

Max user weight 150kg
405mm x 405mm
50mm thick



Yes

Cushion Basic Foam 460mm x
405mm Black Cover

Max user weight 150kg
460mm x 405mm
50mm thick



Yes

Cushion Basic Foam 500mm x
450mm Black Cover

Max user weight 150kg
500mm x 450mm
50mm thick



Yes

iCare Basic Foam Cushion 400
x 400 x 50

Range of sizes and thicknesses

400 x 400 x 50

- » Zip off washable covers
- » Range of sizes available
- » Pressure and heat sensitive
ActiveX™ core

The tandard seat cushion range
features high density visco cores which
respond to pressure in a different way.

Instead of pushing back like standard
foam cushions the Icare cushions will
mould around the shape of your body
taking away points of pressure.



Yes

iCare Basic Foam Cushion 400
x 400 x 100

Range of sizes and thicknesses

400 x 400 x 100

- » Zip off washable covers
- » Range of sizes available
- » Pressure and heat sensitive
ActiveX™ core

The tandard seat cushion range
features high density visco cores which
respond to pressure in a different way.

Instead of pushing back like standard
foam cushions the Icare cushions will
mould around the shape of your body
taking away points of pressure.



Yes

iCare Basic Foam Cushion 450
x 450 x 50

Range of sizes and thicknesses

450 x 450 x 50

- » Zip off washable covers
- » Range of sizes available
- » Pressure and heat sensitive
ActiveX™ core

The tandard seat cushion range
features high density visco cores which
respond to pressure in a different way.

Instead of pushing back like standard
foam cushions the lcare cushions will
mould around the shape of your body
taking away points of pressure.



Yes

iCare Basic Foam Cushion 450
x 450 x 50

Range of sizes and thicknesses

450 x 450 x 100

- » Zip off washable covers
- » Range of sizes available
- » Pressure and heat sensitive
ActiveX™ core

The tandard seat cushion range
features high density visco cores which
respond to pressure in a different way.

Instead of pushing back like standard
foam cushions the Icare cushions will
mould around the shape of your body
taking away points of pressure.



Yes

iCare Basic Foam Cushion 450
x 500 x 50

Range of sizes and thicknesses

450 x 500 x 50

- » Zip off washable covers
- » Range of sizes available
- » Pressure and heat sensitive
ActiveX™ core

The tandard seat cushion range
features high density visco cores which
respond to pressure in a different way.

Instead of pushing back like standard
foam cushions the Icare cushions will
mould around the shape of your body
taking away points of pressure.



Yes

iCare Basic Foam Cushion 450 x 500 x 100

Range of sizes and thicknesses

450 x 500 x 100

- » Zip off washable covers
- » Range of sizes available
- » Pressure and heat sensitive ActiveX™ core

The standard seat cushion range features high density visco cores which respond to pressure in a different way.

Instead of pushing back like standard foam cushions the Icare cushions will mould around the shape of your body taking away points of pressure.



No

Basic Foam Cushion with Gel Insert. (450 x 400 x 50mm)

This foam pressure care cushion is suited to people at low to medium risk of pressure injury that is, damage to the skin and underlying tissue as a result of prolonged sitting without moving. Helps relieve chronic pain and pressure for users up to 110kg. Made from a patented gel bladder encased in resilient PU comfort foam cushion for superior weight equalisation and pressure distribution.

Max user weight: 110Kg
450x400x50mm



No

Basic Foam Cushion with Gel Insert. (450 x 450 x 50mm)

This foam pressure care cushion is suited to people at low to medium risk of pressure injury that is, damage to the skin and underlying tissue as a result of prolonged sitting without moving. Helps relieve chronic pain and pressure for users up to 110kg. Made from a patented gel bladder encased in resilient PU comfort foam cushion for superior weight equalisation and pressure distribution.

Max user weight: 110Kg
450x450x50mm



No

Basic Foam Cushion with Gel Insert. (500 x 450 x 50mm)





This foam pressure care cushion is suited to people at low to medium risk of pressure injury that is, damage to the skin and underlying tissue as a result of prolonged sitting without moving. Helps relieve chronic pain and pressure for users up to 110kg. Made from a patented gel bladder encased in resilient PU comfort foam cushion for superior weight equalisation and pressure distribution.

Max user weight: 110Kg
500x450x50mm



Yes	Contoured Wheelchair Seat Cushion (18" x 16")	An anatomically moulded wheelchair cushion that is used to relieve and prevent back pain, numbness and other discomfort that can be experienced after prolonged sitting in a slung seat wheelchair.	18" x 16" 460 x 410 mm
	Comfort Products		
Yes	Sheepskin Medical Grade	<p>Genuine Australian lambswool overlay Provides additional comfort to existing mattress or chair Soft pile allows air circulation to keep skin dry Natural fibre reduces the risk of skin irritation Machine washable Dimensions are calculated on an average</p>	640mm wide x 960mm Long
Yes	Sheepskin Slippers Small - Pair	<p>Low profile medical sheepskin slipper pressure-redistributing and microclimate controlling properties non-slip suede soles slim and lightweight washable to 80°C</p>	



Yes	Sheepskin Slippers Medium - Pair	Low profile medical sheepskin slipper pressure-redistributing and microclimate controlling properties non-slip suede soles slim and lightweight washable to 80°C	
Yes	Sheepskin Slippers Large - Pair	Low profile medical sheepskin slipper pressure-redistributing and microclimate controlling properties non-slip suede soles slim and lightweight washable to 80°C	
Yes	Sheepskin Heel Protectors	Genuine medical sheepskin Soft pile allows air circulation to keep skin dry Natural fibres reduces the risk of skin irritation Adjustable velcro securing straps Machine washable One size fits all Supplied as a pair.	
Yes	Heel Protector - Bed Heel Elevator	Wedge-shaped contoured foam provides excellent support for the lower leg Rigid foam aids blood flow and reduces pressure	
Static Mattresses			

No

Mattress - Area Care Plus -
Static Foam

The AreaCare Plus is available in Standard and King Single size. The three layered, zoned design of this very high risk pressure surface is ideal for patients with poor mobility or uncertain skin sensitivity.

Three-layer foam pressure surface offers maximum pressure reduction in a three-layer foam surface for higher-risk patients.

A top layer of high density, temperature-sensitive, memory foam conforms to the body's shape, working in combination with the higher resistance support layer to disperse pressure evenly. The top layer also incorporates a special immersion, slow recovery foam foot zone for the ulcer prone heel area. The core layer and edges provide form, strength, and stability. Standard model suitable for users up to 220 kilograms; King Single up to 250 kilograms.

Single (1960 x 865 x
150mm)



No

Mattress - Area Care Pro Static Foam

The AreaCare Pro multi-density foam pressure surface offers maximum pressure reduction in a three layer foam surface for higher risk patients. The three layered, zoned design of this very high risk pressure surface is ideal for heavier patients with poor mobility or uncertain skin sensitivity. Made from premium materials, each foam layer is made of varying density to maximise the redistribution of pressure. The top layer memory foam also incorporates a special immersion, slow recovery foam foot zone for the ulcer prone heel area, improves ventilation to the skins surface and increases the overall weight capacity to 220 kilograms.

Single (1960 x 865 x 150mm)



No

iCare Long Single Mattress - SOFT

For clients who like a soft mattress, or who spend a long time in bed. The feel of the IC25 soft Mattress is best suited to aged or petite patients.

The IC25 Mattresses have almost no partner disturbance or roll together.

Long Single (2030 x 900 x 225)



No

iCare Long Single Mattress -
MEDIUM

Middle of the range, the IC20 proves to be popular for most. The medium feel provides a soft feel but excellent support.

How pressure care works with the IC20: This IC20 medium mattress is heat and pressure sensitive which means that as it gets warmer it becomes softer. So when there is a pressure point developing which usually gets hot and enflamed, this area of the mattress will become softer which allows the blood to flow freely again.

Long Single (2030 x 900
x 200)



No

iCare Long Single Mattress -
FIRM

The IC15 firm mattress offers support and comfort. Firm feel makes it easier for turning.

Long Single (2030 x 900
x 150)



No

iCare King Single Mattress -
SOFT

For clients who like a soft mattress, or who spend a long time in bed. The feel of the IC25 soft Mattress is best suited to aged or petite patients.

King Single (2030 x 1070 x 225)

The IC25 Mattresses have almost no partner disturbance or roll together.



No

iCare King Single Mattress -
MEDIUM

Middle of the range, the IC20 proves to be popular for most. The medium feel provides a soft feel but excellent support.

King Single (2030 x 1070 x 200)

How pressure care works with the IC20: This IC20 medium mattress is heat and pressure sensitive which means that as it gets warmer it becomes softer. So when there is a pressure point developing which usually gets hot and enflamed, this area of the mattress will become softer which allows the blood to flow freely again.



No

No

No

iCare King Single Mattress - FIRM

The IC15 firm mattress offers support and comfort. Firm feel makes it easier for turning.

King Single (2030 x 1070 x 150)



Alternating Pressure Reduction Mattresses

ProCair Mattress Replacement System

ProCair is designed for individuals at the high-risk pressure injury category. Built for optimal patient therapy and comfort, ProCair is now available in Standard and King Single sizes with sealed base to support your infection control strategy.

Weight Capacity (kg)
40-200



ProCair Plus Mattress Replacement System

The ProCair Plus is our premium solution for high-risk patients optimal patient therapy and comfort; fully automated operational features; advanced safety features including back-up battery; and advanced error and service indicators. Available in Standard and King Single models.

King Single (2000 x 1050 x 200mm)
Single (2000 x 880 x 200mm)
Weight Capacity 40-200Kg



No

Pressure Reduction Overlays

Procair Mattress Overlay
System

Functional overlay for pressure prevention
Automatic pressure setting
Used over the existing mattress for basic support
Easy to operate and maintain
40-150 kg
Ideally suited to home care or aged care environments, this product can be used for people at risk to high risk of pressure injury.
Automatic weight sensing technology that delivers fail-safe patient comfort levels.
Conform to all relevant Australian standards.

2000 x 880 x 130 mm
Weight Capacity (kg) 40-150kg



CHSP EQUIP
IN SCOPE LIST

PRODUCT

PRODUCT DESCRIPTION

SPECIFICATIONS

IMAGE

TRANSFER AIDS

Handy Bars

Yes

Transfer Handy Bar

Provides a safe and supportive handhold for transfers in and out of the car
Multifunctional device designed to help users get in and out of a car easier.

Maximum User Weight
Capacity: 320 kg



Leg Lifters

Yes

Leg Lifter (Single) Pelican Brand

Rigid device with loop for foot at one end and two straps at the other for assisting independence to lift a leg into a car or bed.

Overall length: 102 cm



Yes






Dual Handle Leg Lifter




A single leg lifter with 2 handles for users who cannot lift their leg with one hand.

Length 86.5 cm
Weight 110g

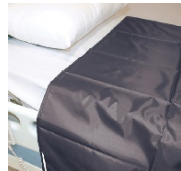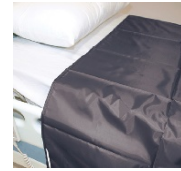

Slide Sheets

Yes	Slide Sheet SMALL - 1000mm x 1450mm	Slippery surface both sides, use folded onto itself to provide low resistance under clients requiring assistance with repositioning. Size: 1m x 1.4m sq for bed transfers	1m x 1.4m vary	Colour may vary	
Yes	Slide Sheet LARGE - 2000mm x 1450mm	Slippery surface both sides, use folded onto itself to provide low resistance under clients requiring assistance with repositioning. Size: 2m x 1.4m sq for bed transfers	2m x 1.4m vary	Colour may vary	
Slide Boards					
Yes	Car Transfer Pad - Pelican	Car slide transfer pad has 2 slip surfaces that decreases friction allowing the client to slide and pivot into the car seat.	Top Pad size 58cm x 58cm		
Yes	The Slider	"Assists patient transfers between wheelchair and car, chair, bed etc Strong perspex construction Rear wheel cutout, covered with a flexible plastic guard, minimises skin damage while transferring over the arch of the wheel Not to be used as a lifting device"	Length 640mm Width 350mm Weight 1.5kg Max User Weight 120kg		
Yes	PVC Angled Transfer Board	A smaller transfer board angled at one end to allow better positioning around the wheel of a wheelchair for side transfers. 15833 Angled Transfer Board (P.V.C.) flexible construction for slight height variation.	SWL 110kg 600mm x250mm x 6mm		

Yes	Curved Transfer Board PVC	PVC transfer board	Weight: 1.95kg. Max. user weight: 152kg	
Turntables				
Yes	Swivel Cushion Rota	<p>This cushioned seat glides you in and out of cars smoothly without struggling. Placed on any seat it will swivel in both directions and may help those who suffer from arthritis, leg or back pain. It's lightweight, portable and consists of a plastic base with a fleece covered cushion. A strap secures the position of the swivel seat.</p>	<p>Diameter: 40cm Height: 70mm Weight: 1.6kg Maximum user weight: 190kg</p>	
Yes	Patient Swivel Seat	<p>Reduces the risk of strain or back injury for the carer and user Assists the user in transfers from a seated position Ideal for transfers in a variety of situations such as car, wheelchair and household areas Swivels 360° Lightweight and portable Smooth and safe turning action Covered in soft fabric upholstery Nonslip underside</p>	<p>Diameter: 39cm Height: 50mm Weight: 0.9kg Maximum user weight: 100kg</p>	

Yes

Handy Soft Turn Pad

Reduces the risk of strain or back injury for the carer and user
Assists the user in transfers from a seated position
Ideal for transfers in a variety of situations such as car, wheelchair and household areas
Swivels 360°
Lightweight and portable
Cotton upholstered seating surface
Nylon swivel mechanism provides smooth and safe turning action
Nonslip underside

Diameter: 45cm

Height: 40mm

Weight: 0.35kg

Maximum user weight:
125kg



CHSP EQUIP
IN SCOPE LIST

PRODUCT

PRODUCT DESCRIPTION

SPECIFICATIONS

IMAGE

DAILY LIVING AIDS

Crockery

Yes

Feeding Cup with Long Spout
and HandlesHomecraft Feeding Cup, Standard
8mm Spout, 250ml, 2/pack

Yes

Plate Incurve Surround
Homecraft

This flexible guard can be fitted to ordinary plates to assist with one-handed eating. Three clips attach it to the rim and create a high inward sloping face to help with food collection and prevent spillage.

Plate size 190 to 254mm external diameter. Dishwasher safe. Melamine material dishwasher safe to 180 degrees F (82°C).



Yes

Plate Manoy Large

Manoy Contoured Plate, Large, 279 x
197mm, White


Cutlery




Yes



SURE GRIP - Teaspoon

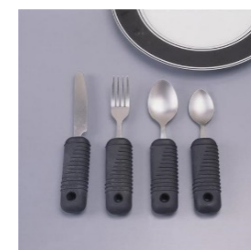

Comfortable ribbed built-up handles help people with weak grasp to hold utensils securely. Stainless steel with rubber handles. Institutional dishwasher safe up to 125 F (52°C)

SURE GRIP Non-bendable Cutlery Dishwasher safe.



Yes	SURE GRIP - Tablespoon	Comfortable ribbed built-up handles help people with weak grasp to hold utensils securely. Stainless steel with rubber handles. Institutional dishwasher safe up to 125 F (52°C)	SURE GRIP Non-bendable Cutlery Dishwasher safe.	
Yes	SURE GRIP - Fork	Comfortable ribbed built-up handles help people with weak grasp to hold utensils securely. Stainless steel with rubber handles. Institutional dishwasher safe up to 125 F (52°C)	SURE GRIP Non-bendable Cutlery Dishwasher safe.	
Yes	SURE GRIP - Knife	Comfortable ribbed built-up handles help people with weak grasp to hold utensils securely. Stainless steel with rubber handles. Institutional dishwasher safe up to 125 F (52°C)	SURE GRIP Non-bendable Cutlery Dishwasher safe.	

Yes	SURE GRIP BENDABLE - Teaspoon	SURE GRIP Cutlery Bendable Utensils are made of stainless steel featuring a 35mm diameter built up handle.	SureGrips Cutlery Bendable Utensils are made of stainless steel featuring a 35mm diameter built up handle. The forks and spoons have a special twist built into the metal shaft that allows them to be bent in any angle for either left or right handed use. Handle length 165mm. Dishwasher safe.	
Yes	SURE GRIP BENDABLE - Tablespoon	SURE GRIP Cutlery Bendable Utensils are made of stainless steel featuring a 35mm diameter built up handle.	SURE GRIP Cutlery Bendable Utensils are made of stainless steel featuring a 35mm diameter built up handle. The forks and spoons have a special twist built into the metal shaft that allows them to be bent in any angle for either left or right handed use. Handle length 165mm. Dishwasher safe.	

Yes	SURE GRIP BENDABLE - Fork	SURE GRIP Cutlery Bendable Utensils are made of stainless steel featuring a 35mm diameter built up handle.	SURE GRIP Cutlery Bendable Utensils are made of stainless steel featuring a 35mm diameter built up handle. The forks and spoons have a special twist built into the metal shaft that allows them to be bent in any angle for either left or right handed use. Handle length 165mm. Dishwasher safe.	
Yes	SURE GRIP BENDABLE - Knife	SURE GRIP Cutlery Bendable Utensils are made of stainless steel featuring a 35mm diameter built up handle. Please note: The knife is not bendable but features a unique rocking action.	SURE GRIP Cutlery Bendable Utensils are made of stainless steel featuring a 35mm diameter built up handle. The forks and spoons have a special twist built into the metal shaft that allows them to be bent in any angle for either left or right handed use. Handle length 165mm. Dishwasher safe.	

Yes

GOOD GRIPS BENDABLE - Teaspoon

Good Grips Bendable Utensils offer forks and spoons that can be adjusted to different angles. The shaft of the utensils features a special twisting mechanism so that users can customize the reach of their forks and spoons. This allows for more comfortable use that puts less stress on the hand. Note: Rocker Knife is not bendable.

6.5 inch stainless steel bendable utensils that feature 1-3/8 inch built up rubber handle
Latex-free, flexible handle features ribbing for easy and adjustable gripping
Ideal for individuals that have a limited range of motion, especially with hand to mouth movements
Available as Fork, Teaspoon, Tablespoon, and Rocker Knife
Utensils are dishwasher safe
Handle length 165mm



Yes

GOOD GRIPS BENDABLE - Tablespoon

Good Grips Bendable Utensils offer forks and spoons that can be adjusted to different angles. The shaft of the utensils features a special twisting mechanism so that users can customize the reach of their forks and spoons. This allows for more comfortable use that puts less stress on the hand. Note: Rocker Knife is not bendable.

6.5 inch stainless steel bendable utensils that feature 1-3/8 inch built up rubber handle
Latex-free, flexible handle features ribbing for easy and adjustable gripping
Ideal for individuals that have a limited range of motion, especially with hand to mouth movements
Available as Fork, Teaspoon, Tablespoon, and Rocker Knife
Utensils are dishwasher safe
Handle length 165mm



Yes

GOOD GRIPS BENDABLE - Fork

Good Grips Bendable Utensils offer forks and spoons that can be adjusted to different angles. The shaft of the utensils features a special twisting mechanism so that users can customize the reach of their forks and spoons. This allows for more comfortable use that puts less stress on the hand. Note: Rocker Knife is not bendable.

6.5 inch stainless steel bendable utensils that feature 1-3/8 inch built up rubber handle
Latex-free, flexible handle features ribbing for easy and adjustable gripping
Ideal for individuals that have a limited range of motion, especially with hand to mouth movements
Available as Fork, Teaspoon, Tablespoon, and Rocker Knife
Utensils are dishwasher safe
Handle length 165mm



Yes

GOOD GRIPS NON
BENDABLE - Knife

Good Grips Bendable Utensils offer forks and spoons that can be adjusted to different angles. The shaft of the utensils features a special twisting mechanism so that users can customize the reach of their forks and spoons. This allows for more comfortable use that puts less stress on the hand. Note: Rocker Knife is not bendable.

6.5 inch stainless steel bendable utensils that feature 1- $\frac{3}{8}$ inch built up rubber handle
Latex-free, flexible handle features ribbing for easy and adjustable gripping
Ideal for individuals that have a limited range of motion, especially with hand to mouth movements
Available as Fork, Teaspoon, Tablespoon, and Rocker Knife
Utensils are dishwasher safe
Handle length 165mm



Yes

GOOD GRIPS WEIGHTED and
BENDABLE - Teaspoon





Good Grips Weighted & Bendable Utensils have 170g of added weight in the built-up handle. This provides more control so that the utensil reaches the mouth easier. The utensil head is bendable. The utensils are ideal for persons with limited hand control, Parkinsons disease, or spasticity.





Available as Fork, Teaspoon, Tablespoon, Souper Spoon, and Rocker Knife
Built-up handle
170g
Institutional dishwasher safe







Yes	GOOD GRIPS WEIGHTED and BENDABLE - Tablespoon	Good Grips Weighted & Bendable Utensils have 170g of added weight in the built-up handle. This provides more control so that the utensil reaches the mouth easier. The utensil head is bendable. The utensils are ideal for persons with limited hand control, Parkinsons disease, or spasticity.	Available as Fork, Teaspoon, Tablespoon, Souper Spoon, and Rocker Knife Built-up handle 170g Institutional dishwasher safe
Yes	GOOD GRIPS WEIGHTED and BENDABLE -Souper spoon	Good Grips Weighted & Bendable Utensils have 170g of added weight in the built-up handle. This provides more control so that the utensil reaches the mouth easier. The utensil head is bendable. The utensils are ideal for persons with limited hand control, Parkinsons disease, or spasticity.	Available as Fork, Teaspoon, Tablespoon, Souper Spoon, and Rocker Knife Built-up handle 170g Institutional dishwasher safe
Yes	GOOD GRIPS WEIGHTED and BENDABLE -Fork	Good Grips Weighted & Bendable Utensils have 170g of added weight in the built-up handle. This provides more control so that the utensil reaches the mouth easier. The utensil head is bendable. The utensils are ideal for persons with limited hand control, Parkinsons disease, or spasticity.	Available as Fork, Teaspoon, Tablespoon, Souper Spoon, and Rocker Knife Built-up handle 170g Institutional dishwasher safe



Yes	GOODGRIPS Knife - Weighted	<p>Good Grips Cutlery Bendable Utensils are made of stainless steel featuring a 35mm diameter built up handle. The Weighted Utensils range are available with 170g of extra weight in the built up handle to provide more control for people with a tremor or limited hand control, Parkinson's disease or spasticity. The knife is weighted but not bendable.</p>	<p>Available as Fork, Teaspoon, Tablespoon, Souper Spoon, and Rocker Knife Built-up handle 170g Institutional dishwasher safe</p>	
	Dressing Aids			
Yes	Button Hook Modular Kings Complete	<p>Assists one-handed dressing Stainless steel wire passes through button hole, gripping the button. It is pulled through with a slight twisting motion Contoured handle</p>	<p>Length 230mm Weight 70g</p>	
Yes	Dressing Stick	<p>Wooden stick with epoxy coated hook at one end and rubber tip at the other. Assists in applying and removing clothing. Can also be used to push or pull objects.</p>	<p>66cm length 5/8"-diameter wooden dowel. Finished with lacquer with reinforced hook.</p>	
Yes	Shoe Horn Long Handled Plastic	<p>Plastic handle designed to assist in putting on shoes especially for people who have difficulty bending and reaching.</p>	<p>Length 60cm. Blue Colour. Includes wrist strap.</p>	

Yes	Stocking Aid	<p>Stocking aid allows the user to don socks and stockings independently. Side notches hold the sock securely, while the strap ends help pull stocking into position efficiently. Built with a flexible plastic gutter for improved maneuverability.</p>	<p>Available in contoured or notched. Straps 36". Gutter length 9".</p>	
Yes	Sock Aid - Single Plastic	<p>Simple aid designed to assist with putting on a stocking especially if bending and reaching the foot is difficult.</p> <p>Cut out side notches designed to hold stocking</p> <p>1 metre length pull straps</p>		
Yes	Stocking Compression Applicator Ezy-As Small	<p>The unique design of the Ezy-As compression garment applicator makes it possible to assist in the fitting of a wide range of compression garments.</p>	Small 28 - 35 cm	
Yes	Stocking Compression Applicator Ezy-As Medium	<p>The unique design of the Ezy-As compression garment applicator makes it possible to assist in the fitting of a wide range of compression garments.</p>	Medium 32 - 40cm	

Yes	Stocking Compression Applicator Ezy-As Large	The unique design of the Ezy-As compression garment applicator makes it possible to assist in the fitting of a wide range of compression garments.	Large 38cm - 47cm	
Yes	Compression Stocking / Garment Applicator - Small		Small	
Yes	Sock Aid - Terry Towelling	Simple aid designed to assist with putting on a sock especially if bending and reaching the foot is difficult. Terry towelling outer fabric for greater hold with smooth nylon inner surface to reduce friction.	Nylon pull straps - 700mm	
<hr/>				
Yes	Household Aids Tap Turner - Crystal	This Crystal Tap Turner will assist you to use the tap on any tap that does not have an 'X' style knob. The Turner clamps onto the tap and then you can use the handle to control the tap. This is a fantastic tool for those who struggle to grasp a tap and close their hand. Instead, this tool gives you a larger area handle to grip on to.	Length 190mm. Weight 112g	



Yes




Tap Turner - Red/Blue



The Homecraft Tapturn Turners are strong moulded turners which simply hook over the top of a 'X' type tap, to make it easy to turn. They are supplied as a red and blue pair for hot and cold taps.

Length 135mm
Width 47mm



Yes	Non Slip Matting	<p>A roll of versatile soft material that can be used on any surface to help stop items from moving, or to protect against scratches.</p> <p>This non-slip matting has a slightly spongy texture which prevents items moving or slipping on almost any surface.</p> <p>Place a piece of matting on your table or bench, under trays or when preparing foods - you can even use it under your floor mat to help stop it sliding. Also great to prevent unwanted scratches on any surface. Try wrapping some around a jar lids to improve your grip when opening. The unique plastic material won't fray and and simply wipes clean with a damp cloth - or throw it in the washing machine.</p>	<p>Comes in a roll (size 1000 mm long by 300 mm wide) and can be cut to size with household scissors.</p>	
Yes	<p>Tray Mobile</p> <p>Trolley Walker - Adjustable Height</p>	<p>Includes clip-on plastic trays that can be removed easily for cleaning.</p> <p>Can be used with or without trays or it can be used as a table.</p> <p>Rails on tray edges prevent items from slipping off.</p> <p>Lockable hand brakes mean this trolley can also be used as a walking aid and</p>	<p>Overall Width: 585mm Overall Depth: 560mm Handle Height Range: 755mm - 905mm Wheel Diameter: 150mm (6") Maximum User Weight: 150kg</p>	

Yes	Traymobile (no brakes)	<p>Ideal for indoor use</p> <p>Two removable trays to assist transporting smaller household items</p> <p>Height adjustable</p> <p>Four swivel castors provide greater manoeuvrability</p>	<p>Frame Height 800-950mm</p> <p>Frame Width 490mm</p> <p>Frame Depth 540mm</p> <p>Max User Weight 125kg</p>	
Yes	<p>Kitchen Aids</p> <p>Plastic Spreading Board</p>	no info available		
Yes	Bread Board	<p>Ideal for people with tremors or impaired sight, the clever design of this cutting aid makes kitchen work much easier. Food is firmly held in the device while a series of rubber stopper feet keep the board from moving around on your kitchen bench or table. Great for bread, cheese, vegetables or other food items, your knife slides between the two guides, making it easy to cut straight and safe, with no risk of injury to fingers or hands. The board can be adjusted for a range of different cutting thickness</p>		

Yes	ETAC Food Prep System	<p>This durable cutting board is ideal for those who have difficulty in gripping, or only have the use of one hand.</p> <p>The clamp can fix items such as tins or bowls, and the stainless steel spikes keep vegetables still when peeling or cutting.</p> <p>It has four non-slip rubber feet to prevent the board from moving.</p> <p>Suction pads are also provided for use on wet surfaces.</p>	<p>Size: 315mm × 295mm.</p> <p>Weight: 1.66kg.</p>	
Yes	Cordless Kettle Tipper	<p>Suits most cordless kettles.</p> <p>There is a metal plate that secures the base of the cordless kettle and a frame with a Velcro strap to keep the upper part of the kettle in place.</p> <p>A cut out at the foot of the tipper allows cups to be positioned as close to the kettle as required.</p> <p>The frame has a built-in back stop to prevent spillage if the kettle is accidentally released during use.</p>	<p>Weight: 1kg</p>	

Yes

Uccello Kettle Tipper

The Uccello Tipping Kettle was designed using the innovative and unique PowerPour technology which means the user can safely pour boiling water without removing the kettle from its cradle by just tilting the kettle.

1.5L capacity

Its ergonomic handle makes it comfortable to use for people with arthritis, grip weakness, or hand injury. The kettle features a wide spout and large hinged lid to easily fill with water. The water level indicator assists the user to ensure it isn't overfilled. The base is non-slip and weighted.

The Kettle has a stainless steel heating element, power-on light and auto shut-off with overheating protection.

This product is available in Black & White or Red & White.



Yes

Jar & Bottle Opener - Mighty Lever

The Jar and Bottle Opener has a flexible stainless steel band which is fitted around a lid, and then clamped tight using the screw at the end of the handle. The handle can then be used

Suitable for tops 20-110mm diameter. Weight 88g.



Yes

Jar Twister

A cone shaped rubber moulding with fluted finger grips on the outside, and ribbing on the inside to give a good grip for opening twist-off jar lids.

Suitable for tops 20 to 85mm diameter.
Weight: 100g



Yes




Zyliss 5 in 1 Can Opener




This versatile opener is the ideal gadget for opening pull tabs on cans, plastic and metal twist-off caps, metal bottle caps, pop-tops and jar lids. No more struggling to get the caps and lids off bottles and jars, this handy gadget has all your needs covered. The soft touch grip makes it comfortable and simple to use.

Dimensions:
12Lx15Wx4H



Personal Care Aids

Yes	Long Handled Toe Wiper	Helps people who have limited access to their feet toes to wash and dry between their toes.	Long bendable wire stem with a plastic built-up easy-grip handle Flat wire frame covered in a white towelling pad for washing between the toes Includes 2 pads Length: 710mm	
Yes	Reaching Aid Economy - 67cm	Trigger style mechanism pick up stick Slip resistant jaw provides a secure hold Hook trigger Magnetic tip	67cm	
Yes	Reaching Aid Economy - 82cm	Trigger style mechanism pick up stick Slip resistant jaw provides a secure hold Hook trigger Magnetic tip	82cm	

Yes	Reaching Aid Squeeze - short	<p>Slip resistant rotating jaw designed to pick up objects at different angles</p> <p>Simple hand trigger operation</p> <p>Magnetic tip</p> <p>Wheelchair/walking frame clip supplied.</p>	10204/65 (65 cm)	
Yes	Reaching Aid Squeeze - long	<p>Slip resistant rotating jaw designed to pick up objects at different angles</p> <p>Simple hand trigger operation</p> <p>Magnetic tip</p> <p>Wheelchair/walking frame clip supplied.</p>	10204/81 (81 cm)	
Yes	Bottom Wiper	<p>The Homecraft Toilet Aid provides assistance with independent hygiene.</p> <p>The aid is designed with a curved metal handle and a molded-plastic head that gives the user the proper reach. The handle's recessed serrations allow it to securely grip paper.</p>	Length 25.4cm. Weight 71g.	

Yes

Long Handled Sponge

As you get older, it can be difficult to bend down to wash your feet or to lift them up so you can wash them. Doing so might even put you at risk for a slip and fall in the shower. This long-handled bath sponge has a wire that is strong, but also easily bendable to be rotated, twisted, and moved in any position that helps you accomplish optimal cleaning without bending! It can be used by seniors and people with disabilities to bathe independently, by caretakers to easily clean residents, or for household cleaning tasks.

Length 24". Sponge head 3.75" x 5".



CHSP EQUIP
IN SCOPE LIST

PRODUCT

PRODUCT DESCRIPTION

SPECIFICATIONS

IMAGE

LIFTERS AND HOISTS

Lifters

No

Patient Lifter - 180kg

The LiftAid 180 is a lightweight, easy to operate mobile lifter that allows patients to be lifted and transferred from floor, bed or chair with ease and safety.

- 3 height settings - tool-less adjustment
- Quick release hook system
- Electrically operated base
- Handheld controller
- Emergency stop
- Manual emergency lowering
- Optional Pivot Frame
- Compatible with range of slings

- Width: 61.5cm
- Height: 43cm - 181cm (3 height settings)
- Length: 109cm
- Product Weight: 35.3kg
- Warranty: Lifter & Actuator - 2 years | Battery - 1 year
- S.W.L: 180kg
- Max Spread: 101cm
- Underbed Clearance: 12.5cm



No

Patient Lifter - 320kg

The LiftAid 320 is a heavy-duty mobile lifter ideal for the majority of common lifting situations.

Patients can be lifted and transferred from floor, bed or chair. Fold-down armrests assist with gait training during rehabilitation or for general toileting.

Adjustable height frame allows overall mast height be set to match environmental and lifting needs.

Features the advanced Linak® Jumbo Care Diagnostic Control System for intelligent self-diagnostic system capabilities. Low lifting height and soft start/stop function,

3 height settings with tool-less adjustment

quick release hook system for easy sling bar installation

swivel-style front and rear castors; foot-operated brake locks on rear

3 height settings with tool-less adjustment

Quick release hook system for easy sling bar installation

Swivel style front and rear castors; foot operated brake locks on rear

Electrically operated 'U' shaped base with wide flanged front legs for enhanced stability and access

Length Overall - 1410 mm

Width Overall 7950 mm - 1190 mm

Height Overall 2180 - mm

Height of Lift 650 mm - 1840 mm

Clearance Under Bed 16 - mm

Diameter of Front Wheel - 105 mm

Diameter of Rear Wheel - 127 mm

Load Capacity - 320 kg



Lifter/Hoist Slings

No

Sling General Purpose - Sml

An easy fitting general purpose sling providing good support for patients with basic upper body control. Suits general lifting from the floor, bed or chair. Versatile enough for general patient transfers, toileting or hygiene care. 3 loop support strap system controls the positioning of the head and leg supports
Suits the majority of lifters on the market
Brightly coloured handles and adjusters secured to sling for easy access and visibility
Quick-drying mesh fabric machine washable up to 85°C for microbiological control
Heavy-duty, 3 mm quilted cushioning foam in leg sections for greater comfort

840 x 1160 x 490 mm (W x L x D)
Capacity - 205 kg



No

Sling General Purpose - Med

An easy fitting general purpose sling providing good support for patients with basic upper body control. Suits general lifting from the floor, bed or chair. Versatile enough for general patient transfers, toileting or hygiene care. 3 loop support strap system controls the positioning of the head and leg supports
Suits the majority of lifters on the market
Brightly coloured handles and adjusters secured to sling for easy access and visibility
Quick-drying mesh fabric machine washable up to 85°C for microbiological control
Heavy-duty, 3 mm quilted cushioning foam in leg sections for greater comfort

890 x 1200 x 545 mm (W x L x D) Capacity - 205kg





No



Sling General Purpose - Lge



An easy fitting general purpose sling providing good support for patients with basic upper body control. Suits general lifting from the floor, bed or chair. Versatile enough for general patient transfers, toileting or hygiene care. 3 loop support strap system controls the positioning of the head and leg supports Suits the majority of lifters on the market Brightly coloured handles and adjusters secured to sling for easy access and visibility Quick-drying mesh fabric machine washable up to 85°C for microbiological control Heavy-duty, 3 mm quilted cushioning foam in leg sections for greater comfort

1135 x 1325 x 640 mm
(W x L x D)
Capacity - 205kg



No	Sling General Purpose - Xlge	<p>An easy fitting general purpose sling providing good support for patients with basic upper body control. Suits general lifting from the floor, bed or chair. Versatile enough for general patient transfers, toileting or hygiene care. 3 loop support strap system controls the positioning of the head and leg supports</p> <p>Suits the majority of lifters on the market</p> <p>Brightly coloured handles and adjusters secured to sling for easy access and visibility</p> <p>Quick-drying mesh fabric machine washable up to 85°C for microbiological control</p> <p>Heavy-duty, 3 mm quilted cushioning foam in leg sections for greater comfort</p>		
	Sling General Purpose Head Support Sml	<p>General Purpose with Head Support is an easy fitting general purpose sling ideally suited for patients with less control. Suits general lifting from the floor, bed or chair. Versatile enough for general patient transfers, toileting or hygiene care. Suitable for users up to 205 kilograms (Small and Medium) or 300 kilograms (Large, Extra Large and Extra, Extra Large).</p> <p>For Yoke frame lifters only</p>	900, 230, 450, 790, 965 mm, 205 kg	

No	Sling General Purpose Head Support Med	General Purpose with Head Support is an easy fitting general purpose sling ideally suited for patients with less control. Suits general lifting from the floor, bed or chair. Versatile enough for general patient transfers, toileting or hygiene care. Suitable for users up to 205 kilograms (Small and Medium) or 300 kilograms (Large, Extra Large and Extra, Extra Large).	1010, 230, 490, 890, 1060 mm, 205 kg	
		For Yoke frame lifters only		
No	Sling General Purpose Head Support Lge	General Purpose with Head Support is an easy fitting general purpose sling ideally suited for patients with less control. Suits general lifting from the floor, bed or chair. Versatile enough for general patient transfers, toileting or hygiene care. Suitable for users up to 205 kilograms (Small and Medium) or 300 kilograms (Large, Extra Large and Extra, Extra Large).	1080, 250, 480, 990, 1180 mm, 300 kg	
		For Yoke frame lifters only		

No	Sling General Purpose Head Support XLge	<p>General Purpose with Head Support is an easy fitting general purpose sling ideally suited for patients with less control. Suits general lifting from the floor, bed or chair. Versatile enough for general patient transfers, toileting or hygiene care. Suitable for users up to 205 kilograms (Small and Medium) or 300 kilograms (Large, Extra Large and Extra, Extra Large).</p> <p>For Yoke frame lifters only</p>	
No	Sling Hygiene Sml	<p>An easy fitting sling specifically designed for toilet transfer with open area from mid back to mid thigh. open area from mid back to mid thigh for hygiene access</p> <p>buckle waist strap for added security and increased infection control.</p> <p>Generous padding under arms</p> <p>3 loop support strap system controls positioning of head and leg supports.</p> <p>Suits majority of lifters on the market.</p> <p>brightly coloured handles and adjusters, secured to sling for easy access and visibility</p> <p>quick drying mesh fabric; machine washable up to 85°C for microbiological control</p> <p>manufactured from the highest quality material.</p>	

No

Sling Hygiene Med

An easy fitting sling specifically designed for toilet transfer with open area from mid back to mid thigh. open area from mid back to mid thigh for hygiene access
buckle waist strap for added security and increased infection control.
Generous padding under arms
3 loop support strap system controls positioning of head and leg supports.
Suits majority of lifters on the market.
brightly coloured handles and adjusters, secured to sling for easy access and visibility
quick drying mesh fabric; machine washable up to 85°C for microbiological control
manufactured from the highest quality material.

Maximum user weight
205kg



No

Sling Hygiene Lge

An easy fitting sling specifically designed for toilet transfer with open area from mid back to mid thigh. open area from mid back to mid thigh for hygiene access
buckle waist strap for added security and increased infection control.
Generous padding under arms
3 loop support strap system controls positioning of head and leg supports.
Suits majority of lifters on the market.
brightly coloured handles and adjusters, secured to sling for easy access and visibility
quick drying mesh fabric; machine washable up to 85°C for microbiological control
manufactured from the highest quality material.

Maximum user weight
300kg



No

Sling Hygiene XLge

An easy fitting sling specifically designed for toilet transfer with open area from mid back to mid thigh. open area from mid back to mid thigh for hygiene access
buckle waist strap for added security and increased infection control.
Generous padding under arms
3 loop support strap system controls positioning of head and leg supports.
Suits majority of lifters on the market.
brightly coloured handles and adjusters, secured to sling for easy access and visibility
quick drying mesh fabric; machine washable up to 85°C for microbiological control
manufactured from the highest quality material.

Maximum user weight
300kg



No

Sling Hygiene With Head Support Sml

Open area from mid back to mid thigh for hygiene access
inbuilt head support for patients with less control
buckle waist strap for added security and increased infection control
generous padding under arms
3 loop support strap system controls positioning of head and leg supports
brightly coloured handles and adjusters, secured to sling for easy access and visibility
quick drying mesh fabric; machine washable up to 85°C for microbiological control
suits majority of lifters on the market
manufactured from the highest quality material, in compliance with industry standards (AS/NZS ISO 10535:2011)

Maximum user weight 205kg



No

Sling Hygiene With Head
Support Med

Open area from mid back to mid thigh
for hygiene access
inbuilt head support for patients with
less control
buckle waist strap for added security
and increased infection control
generous padding under arms
3 loop support strap system controls
positioning of head and leg supports
brightly coloured handles and
adjusters, secured to sling for easy
access and visibility
quick drying mesh fabric; machine
washable up to 85°C for microbiological
control
suits majority of lifters on the market
manufactured from the highest quality
material, in compliance with industry
standards (AS/NZS ISO 10535:2011)

Maximum user weight
205kg





No

Sling Hygiene With Head
Support Lge

Open area from mid back to mid thigh
for hygiene access
inbuilt head support for patients with
less control
buckle waist strap for added security
and increased infection control
generous padding under arms
3 loop support strap system controls
positioning of head and leg supports
brightly coloured handles and
adjusters, secured to sling for easy
access and visibility
quick drying mesh fabric; machine
washable up to 85°C for microbiological
control
suits majority of lifters on the market
manufactured from the highest quality
material, in compliance with industry
standards (AS/NZS ISO 10535:2011)

Maximum user weight
300kg



No	Sling Hygiene With Head Support XLge	<p>Open area from mid back to mid thigh for hygiene access</p> <p>inbuilt head support for patients with less control</p> <p>buckle waist strap for added security and increased infection control</p> <p>generous padding under arms</p> <p>3 loop support strap system controls positioning of head and leg supports</p> <p>brightly coloured handles and adjusters, secured to sling for easy access and visibility</p> <p>quick drying mesh fabric; machine washable up to 85°C for microbiological control</p> <p>suits majority of lifters on the market</p> <p>manufactured from the highest quality material, in compliance with industry standards (AS/NZS ISO 10535:2011)</p>	Maximum user weight 300kg	
	Stand Aid Slings			
No	Sling - Standaid Torso Medium	<p>Fully padded, including soft padding under the arms</p> <p>Various attachments</p> <p>Built-in Keyhole Plate Safety Straps (Applicable for the version of this product which uses Keyhole Plates)</p> <p>A comfortable, well padded, high back Stand Up Sling with a strong wide</p>	<p>Chest Size 85cm - 110cm</p> <p>Load capacity 200kg</p>	

CHSP EQUIP
IN SCOPE LIST

PRODUCT

PRODUCT DESCRIPTION

SPECIFICATIONS

IMAGE

MOBILITY AID ACCESSORIES

IV POLES

Yes

IV Stand Mobile

2 Hook design
560mm nylon base
Chromed steel pole
50mm castors

Adjustable height 1240-
2130mm



ROLLATOR ACCESSORIES

Yes

Rollator Ski Replacement (x2)

Pair of Rollator Skis to replace lost or
damaged skis.

To fit 25mm Rollator
frames



OXYGEN CARRIER

Yes

Oxygen Bottle Carrier For Wheelchair

Oxygen bottle carrier frame to attach to
wheelchair. Contact store for availablity
or suitability to client's chair.



Yes

Oxygen Carrier for 4 Wheel Walker

Wire frame designed to fit under the seat of the Rollator Deluxe range.

Suits oxygen bottles up to 125 mm diameter (C size)

Allows the under seat bag or basket to remain on Rollator

Angled for optimal oxygen supply & security



Yes

STUMP SUPPORTS
Fracture Board

Designed to support an extended leg and is normally used in a wheelchair. The Board can be used for either the left or right leg. Made from a strong plastic that is easy to wipe clean, to help Infection Control. This plastic is virtually maintenance free and will not splinter. The webbing straps can be removed for washing or replacing if necessary.

Length 112cm
Width 40cm



CHSP EQUIP
IN SCOPE LIST

PRODUCT

PRODUCT DESCRIPTION

SPECIFICATIONS

IMAGE

Yes

Exercise Equipment**Exercise Pedals**

Lightweight and portable, the versatile product can be used for a variety of alternative exercises, for both lower and upper body. The compact design makes it easy to move the exerciser from room to room, to take outside in the fresh air or even to take when travelling. The exerciser is designed to be placed on the floor when exercising with the feet, or placed on a table top when exercising with the hands. It features both forwards and backwards movement, to ensure you are working all muscle groups. It also offers an easily adjustable resistance level - simply turn the knob to either increase or decrease the resistance at will.

This product weighs 2.5 kilograms and measures 280 mm high, 445 mm wide and 530 mm deep.



Yes

Balance Pad

Balance Pad ultrasoft speciality foam is both water- and dirt-resistantApplication:Balance training, co-ordination and reaction training, mobility training, standing stability training, functional retraining of the muscles in the lower extremities, motor-skill training to maintain balance, e.g. for age related physiological changes to the muscles and skeletal frame.The Physiomed Balance Pad - several can be laid on top of each other to increase the "wobble effect". The closed-cell robust and hygienic. They are multi-functional in use and are ideally suited for use as swimming floats.

Length 50cm
41cm
0.7kg





Width
Weight



Pain Management

Yes	Tens Unit Spectrum	<p>This easy to use, digital TENS machine offers natural, non-invasive pain relief without the use of chemicals or drugs. Transcutaneous Electrical Nerve Stimulation (TENS) delivers tiny electrical impulses to the nerve fibres that lie underneath the skin's surface, providing pain relief by blocking pain signals to the brain via the spinal cord and peripheral nervous system.</p> <p>Large liquid crystal display and easy operation press buttons, ideal for people who have difficulty reading smaller panel inscriptions or turning smaller dials. The unit offers 32 stimulation modes for pain relief and muscle training. There are 5 adjustable TENS modes, 3 adjustable EMS modes and 24 professional presets for TENS and EMS available.</p>	<ul style="list-style-type: none"> 1 x TENS machine unit 1 x D size 9-volt battery 1 x AC power adaptor 2 x electrode leads 4 x self adhesive electrodes (pack of 4) 1 x instruction manual 1 x hard plastic carrying case
Yes	Weights		
	HAND WEIGHT Plastic COATED .5KG	<p>These dumbbells are made of steel but coated in a soft vinyl (PVC) easy-clean coating. Can be used during work-outs to increase strength training, or used alone as part of a weight training program.</p>	<p>Dumbbell, Vinyl Coated, Purple 0.5kg</p>



Yes	HAND WEIGHT Plastic COATED 1.0KG	These dumbbells are made of steel but coated in a soft vinyl (PVC) easy-clean coating. Can be used during work-outs to increase strength training, or used alone as part of a weight training program.	Dumbell, Vinyl Coated, Pink 1.0 kg	
Yes	HAND WEIGHT Plastic COATED 1.5KG	These dumbbells are made of steel but coated in a soft vinyl (PVC) easy-clean coating. Can be used during work-outs to increase strength training, or used alone as part of a weight training program.	Dumbell, Vinyl Coated, Green 1.5kg	
Yes	HAND WEIGHT Plastic COATED 2.0KG	These dumbbells are made of steel but coated in a soft vinyl (PVC) easy-clean coating. Can be used during work-outs to increase strength training, or used alone as part of a weight training program.	Dumbell, Vinyl Coated, Mauve 2.0kg	
Exercise Media				
Yes	Theraband Yellow 5.5m	TheraBand Resistance Exercise Bands are best used before working out for preventative stretches that prime muscles for strenuous exercises. They are an essential part of any workout or athletic activity.	Yellow thin 5.5m	

Yes

Theraputty Tan Extra Soft 85g

Metron Therapeutic Putty works to strengthen your hands and can help prevent future hand and wrist problems.

It is also perfect for hand grip exercises, stress relief exercises, wrist exercises, and upper arm exercises. It is superior to squeeze balls, hand weights, hand springs, or any other exercise/stress relief product.

Tan Extra Soft 85g



Yes

Theraputty Yellow Soft 85g

Metron Therapeutic Putty works to strengthen your hands and can help prevent future hand and wrist problems.

It is also perfect for hand grip exercises, stress relief exercises, wrist exercises, and upper arm exercises. It is superior to squeeze balls, hand weights, hand springs, or any other exercise/stress relief product.

Yellow Soft 85g



Yes

Theraputty Red Medium 85g

Metron Therapeutic Putty works to strengthen your hands and can help prevent future hand and wrist problems.
It is also perfect for hand grip exercises, stress relief exercises, wrist exercises, and upper arm exercises. It is superior to squeeze balls, hand weights, hand springs, or any other exercise/stress relief product.

Red Medium 85g



Yes

Theraputty Green Firm 85g

Metron Therapeutic Putty works to strengthen your hands and can help prevent future hand and wrist problems.
It is also perfect for hand grip exercises, stress relief exercises, wrist exercises, and upper arm exercises. It is superior to squeeze balls, hand weights, hand springs, or any other exercise/stress relief product.

Green Firm 85g



Yes

Weights

Ankle Weight - Vinyl 0.9kg
(This is a substitute for the 1.0kg ankle weight which is currently not available)

Vinyl Cuff Weights with a long velcro closure strap assures a secure and comfortable fit even during the most strenuous exercise program.

White 0.9kg Ankle Weight with velcro fastening.



CHSP EQUIP
IN SCOPE LIST

PRODUCT

PRODUCT DESCRIPTION

SPECIFICATIONS

IMAGE

BRACES AND SUPPORTS

BRACES

Yes

Knee Brace - Futuro Sport
Adjustable

Walk. Jog. Run. The FUTURO Sport Adjustable Knee Support keeps you moving with compression, strength and support where you need it. It also helps provide protection on the track, court, or green. With this adjustable knee wrap, there is no excuse not to be active anymore.

Special Features:
Adjustable straps for customised fit and support
Use For:
Sprains, strains, arthritis
Product Type:
Wrap Around
Features and Benefits:

Open knee cap design for added support
Side ventilation for breathability and reduced bunching behind knee
Durable neoprene-blend material for strong support, soothing warmth, and comfort
Easy to put on and take off

One Size Fits Most,
Measure knee
circumference: 33cm -
45.7cm



Yes

Wrist Braces Roylan D Ring Small
Left

Wrist Circumference
14.6cm to 16.5cm
Length 16.5cm



The right choice for high-compliance treatment of arthritis, carpal tunnel syndrome, tendonitis and other wrist disorders. Prevents unwanted flexion or extension, yet allows unobstructed motion of the fingers, thumb. Non-elastic circumferential design for easy application and for firm, evenly distributed pressure without restricting circulation. Protective padding lined with cotton stockinette (except for white) D-ring closure system for simple application, removal and customized tension control. A double layer of stockinette under closures cushions straps. Forearm length is proportional to brace size. Removable, adjustable aluminum bar for additional support. Machine wash and air dry. Latex free.

Yes

Wrist Braces Roylan D Ring Small
Right

Wrist Circumference
14.6cm to 16.5cm
Length 16.5cm



The right choice for high-compliance treatment of arthritis, carpal tunnel syndrome, tendonitis and other wrist disorders. Prevents unwanted flexion or extension, yet allows unobstructed motion of the fingers, thumb. Non-elastic circumferential design for easy application and for firm, evenly distributed pressure without restricting circulation. Protective padding lined with cotton stockinette (except for white) D-ring closure system for simple application, removal and customized tension control. A double layer of stockinette under closures cushions straps. Forearm length is proportional to brace size. Removable, adjustable aluminum bar for additional support. Machine wash and air dry. Latex free.

Yes

Wrist Braces Roylan D Ring
Medium Left

Wrist Circumference
17.1cm to 19.1cm Length
17.8cm

The right choice for high-compliance treatment of arthritis, carpal tunnel syndrome, tendonitis and other wrist disorders. Prevents unwanted flexion or extension, yet allows unobstructed motion of the fingers, thumb. Non-elastic circumferential design for easy application and for firm, evenly distributed pressure without restricting circulation. Protective padding lined with cotton stockinette (except for white) D-ring closure system for simple application, removal and customized tension control. A double layer of stockinette under closures cushions straps. Forearm length is proportional to brace size. Removable, adjustable aluminum bar for additional support. Machine wash and air dry. Latex free.



Wrist Braces Roylan D Ring
Medium Right

Wrist Circumference
17.1cm to 19.1cm Length
17.8cm

The right choice for high-compliance treatment of arthritis, carpal tunnel syndrome, tendonitis and other wrist disorders. Prevents unwanted flexion or extension, yet allows unobstructed motion of the fingers, thumb. Non-elastic circumferential design for easy application and for firm, evenly distributed pressure without restricting circulation. Protective padding lined with cotton stockinette (except for white) D-ring closure system for simple application, removal and customized tension control. A double layer of stockinette under closures cushions straps. Forearm length is proportional to brace size. Removable, adjustable aluminum bar for additional support. Machine wash and air dry. Latex free.



Yes

Wrist Braces Roylan D Ring Large
Left

Wrist Circumference
19.7cm to 21.6 cm
Length 19cm

The right choice for high-compliance treatment of arthritis, carpal tunnel syndrome, tendonitis and other wrist disorders. Prevents unwanted flexion or extension, yet allows unobstructed motion of the fingers, thumb. Non-elastic circumferential design for easy application and for firm, evenly distributed pressure without restricting circulation. Protective padding lined with cotton stockinette (except for white) D-ring closure system for simple application, removal and customized tension control. A double layer of stockinette under closures cushions straps. Forearm length is proportional to brace size. Removable, adjustable aluminum bar for additional support. Machine wash and air dry. Latex free.



Yes

Wrist Braces Roylan D Ring Large
Right

Wrist Circumference
19.7cm to 21.6 cm
Length 19cm

The right choice for high-compliance treatment of arthritis, carpal tunnel syndrome, tendonitis and other wrist disorders. Prevents unwanted flexion or extension, yet allows unobstructed motion of the fingers, thumb. Non-elastic circumferential design for easy application and for firm, evenly distributed pressure without restricting circulation. Protective padding lined with cotton stockinette (except for white) D-ring closure system for simple application, removal and customized tension control. A double layer of stockinette under closures cushions straps. Forearm length is proportional to brace size. Removable, adjustable aluminum bar for additional support. Machine wash and air dry. Latex free.



ORTHOSSES

Yes

Rebound Foot Up Ankle Cuff S/M

Rebound Foot-Up is a lightweight ankle foot orthosis that offers dynamic and discreet support for drop foot or related disorders requiring dorsiflexion support

Anke Cuff 18-25cm Foot wrap, strap and bolt 20-28cm

Contoured to the shape of the leg, it is comfortable to wear and easy to use

Rebound Foot-Up provides visible improvement in gait by supporting the foot the moment it is raised



Yes

Rebound Foot Up Ankle Cuff L/XL

Rebound Foot-Up is a lightweight ankle foot orthosis that offers dynamic and discreet support for drop foot or related disorders requiring dorsiflexion support

Anke Cuff 26-33cm Foot wrap, strap and bolt 20-28cm

Contoured to the shape of the leg, it is comfortable to wear and easy to use

Rebound Foot-Up provides visible improvement in gait by supporting the foot the moment it is raised



Yes

Ankle Foot Orthosis Foot-Up
Shoeless Wrap Universal

Rebound Foot-Up Foot Wrap can be used when barefoot or in open shoes or sandals. The adjustable elastic strap of the Rebound Foot Up is attached to the Foot Wrap using a nut and bolt fastened through the eyelets.

The Rebound Foot-Up Ankle Cuff and Rebound Foot-Up Foot Wrap are made from breathable material that can be worn comfortably for long periods of time.

The Rebound Foot-Up is hand-washable.

One size cut to size (20-28cm)



CHSP EQUIP
IN SCOPE LIST

PRODUCT

PRODUCT DESCRIPTION

SPECIFICATIONS

IMAGE

HOME ACCESS


Portable Ramps

Yes

Multipurpose Ramp 870mm

Fibreglass ramp that has a concertina folding action for storage and transport. 875 mm length.

835mm width
875mm length
5.6kg

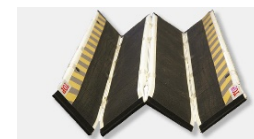


Yes

DecPac Portable Ramp 700mm

Fibreglass ramp that has a concertina folding action for storage and transport. 700mm length.

745mm width
700mm length
4.7kg

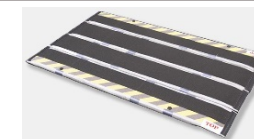


Yes

DecPac Portable Ramp 1200mm

Fibreglass ramp that has a concertina folding action for storage and transport. 1200mm length.

750mm width
1200mm length
8.0kg

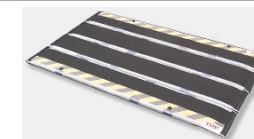


Yes

DecPac Portable Ramp 1200mm (with EDGE)

Fibreglass ramp that has a concertina folding action for storage and transport. 1200mm length with EDGE.

820mm width
1200mm length
9.7kg

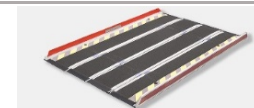


Yes

DecPac Portable Ramp 1650mm

Fibreglass ramp that has a concertina folding action for storage and transport. 1650mm length.

750mm width
1650mm length
10.8kg

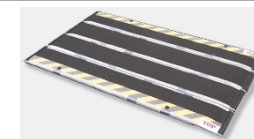



Yes

DecPac Portable Ramp 1650mm (with EDGE)

Fibreglass ramp that has a concertina folding action for storage and transport. 1650mm length with EDGE.

820mm width
1650mm length
13kg

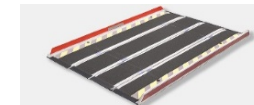

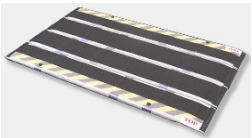Yes	DecPac Portable Ramp 2000mm	Fibreglass ramp that has a concertina folding action for storage and transport. 2000mm length.	750mm width 2000mm length 11.7kg	
Yes	DecPac Portable Ramp 2000mm (with EDGE)	Fibreglass ramp that has a concertina folding action for storage and transport. 2000mm length with EDGE.	820mm width 2000mm length 13.7kg	